Ultrasonic Toothbrush

ProBrush 9.0

medica+

JAPAN TECHNOLOGY

medica+

JAPAN TECHNOLOGY

Instruction Manual Ultrasonic ToothBrush

ProBrush 9.0

Content

Safety Tips	01 - 04
Fitting	05
Operator Panel Instruction	06 - 07
Function Instruction	08 - 10
Tooth-Brushing	11
Timing	12
Install & Use	13 - 14
Charging & Instruction	15 - 16
Clean & Maintenance	17
Environmental Protection	18 - 19
Fault Dignosis	20- 21
Specification	22

Safety tips

🚫 Prohibition!

■ Don't charge in the bathroom! Don't place the charger above or near water container (such as tub pool, etc).

■ Battery can not be replaced by yourself, can only replaced by authorized maintenance center. Manufacturer takes no responsibility for any accidents during replacement by users themselves. See page 13 in case of special reasons to be required to remove by self.

■ In order to avoid danger, once the power cord damages, please stop use this product at once.

⚠ Caution!

■ To ensure safety, 8-12 years old children, mental disorders and disabilities must be supervised or guided when they using this product.

⚠ Attention!

■ It is a normal phenomenon to feel stimulation when using this product at the first time. Please stop using this product and consult your dentist if feeling continuous stimulation in two weeks.

- Please stop using this product and contact your dentist if feeling pain when using it.
- It is normal that the main body is a little warm when charging, please cut off power supply when feeling very hot, and check it in the customer service center.

Safety tips

- All the product accessories should be used our corresponding accessories, like charger, brush head,charging line and adaptor which pass certification required in this country.
- The product is just only designed to clean your teeth, gum and tongue. Please do not use it for any other purposes. This product is for personal product, not for the use of dental clinics or institutions.
- People who just had the oral or periodontal surgery, or are under the treatment within two months,or there are other symptoms(including people whose dentures are week or gums are massively bleeding, etc.) ,or crowds whose oral is sensitive and vulnerable, please stop to use the product and consult your dentist.
- The toothbrush should only be used at room temperature and should not be exposed to the sun. It should be kept in airy place.
- Please cut off the power of charger when not charging .unplug from socket before . Please do not clean toothbrush head, toothbrush holder, charger, charging line with disinfectant in the dishwasher, washing machine and ultrasonic

cleaner. See page 11 for details of cleaning & care.


■ Please don't leak the gelatin, oil, glue or stimulation chemical liquids on the product surface.

Regularly clean power plug, preventing dust and dirt accumulation.

■ The suggested lifespan should be 3 years. Please don't use over 3 years. You should take out the battery form toothbrush before scrapped. If you want to learn how to dismantle battery please refer to Environmental Protection & Recycle at page 13.


■ When the power is on, the charger should not be in direct contact with or near metal objects, such as keys, coins, cell phone cases, metal caps, etc. The metal body can cause abnormality of wireless inductive charger or even damage the charger.

Fitting

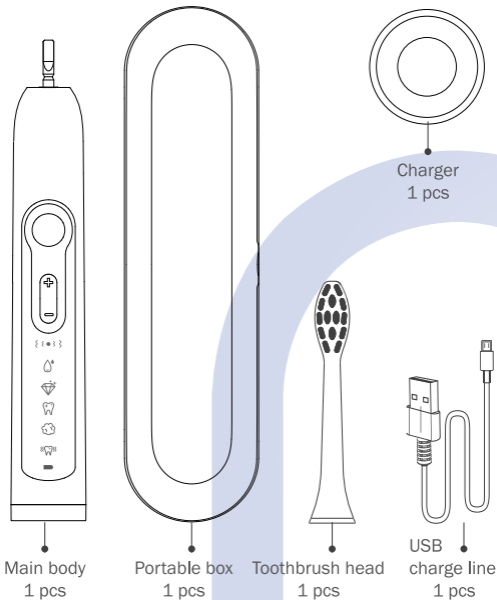* Please carefully check the fitting and its number

Operator panel instruction

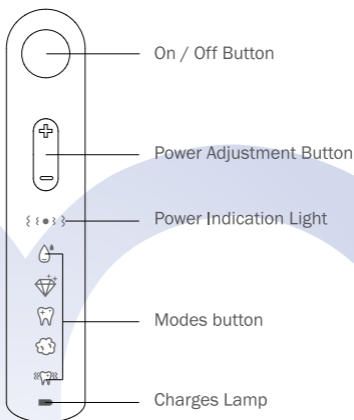■ ON/OFF BUTTON

Short press to start, choose the mode you need. Stop the toothbrush by pressing ON/OFF button for 4 seconds.

■ POWER ADJUSTMENT BUTTON

You can press the "+" or button to change the motor's vibration speed as your own preference to reach the power you need.

■ POWER INDICATION LIGHT

Display current power of toothbrush.
Weak, medium, strong 3 types power display.

■ MODES BUTTON

Display current choosed mode:

-  Clean
-  Polish
-  White
-  Sensitive
-  Massage

■ CHARGING LAMP

Display the condition of charging, full charged, low power.

Function description

■ ON/OFF BUTTON

Short press the On/Off button to turn on the toothbrush and it automatically enters the mode you use last time. The strength defaults to the strong power. Users can adjust the brush power as requirement. Turn off by pressing On/Off 4 seconds.

■ MODES BUTTON

Start the toothbrush according to individual demand to choose the suitable mode, and select mode when the light turning on. The user can switch to other mode at any time during the processing. Short pressing the start key can jump into the next mode. After selection, you can enter the brushing experience.

■ SMART REMINDER

For dental and oral health, the dentist recommends two minutes of brushing time. Every 30 seconds in the middle that your toothbrush will remind you it needs to change direction when there is a pause. (See page 8 for brushing teeth details)

■ INTELLIGENT MEMORY

According to your personal choice, the toothbrush will remember the pattern and the intensity you used last time, then start your brushing experience directly.

■ INTENSITY ADJUSTING

When the toothbrush is turned on, press «—» to weaken the strength. When the indicator light shows {••} is the maximum strength, when the indicator light shows {•} is medium strength, the indicator light shows • is the weakest strength. You also can press «+» to increase strength.

■ MODES

💧 **Clean: 38000 times/min**

High-frequency powerful vibration. This mode is the strongest cleaning mode of this sonic toothbrush. These people are recommended who have used sonic toothbrushes or use sonic toothbrushes frequently.

💎 **Polish: 1000 times/min**

High-frequency weak vibration. This mode is for those who use a sonic toothbrush for the first time or as an adaptive phase for those who just use a sonic toothbrush, or

for these people who don't like too strong power. This mode can help polish your teeth as shiny diamond with high frequency and suitable power.

 **White: 28000-38000 times/min**

High-low frequency rhythm fast vibration to help clean the teeth easily and removes the food debris perfectly.


 **Sensitive: 28000-38000 times/min**

Designed for the sensitive teeth or oral cavity people. Vibration from slow to fast, weak to strong, give the teeth a transitional adaptation process.

 **Massage: 28000-38000 times/min**

The vibration frequency are gradients, from strong to weak, then become strong, cycle like this. It gives good massage to the teeth and gum.

Tooth-brushing

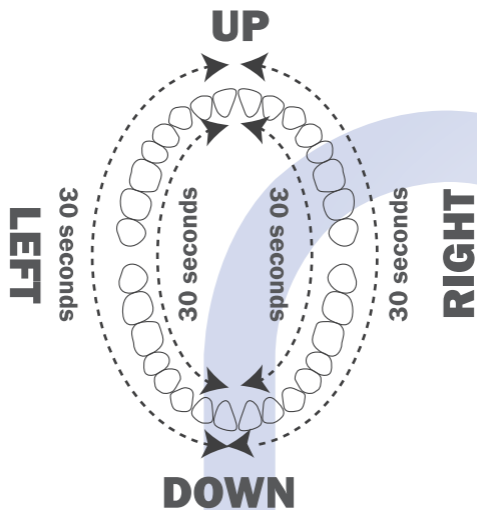
Timing

■ We suggest brushing your teeth twice a day and 2 minutes for each time.


■ You will feel a temporary stop every 30 seconds during use. This aims to reminding you that this part has achieved its best cleaning effect, you need change brush area (Suggest brushing C as showed above)

1. Lateral teeth upper left & lower left
2. Inside teeth upper left & lower left
3. Inside teeth upper right & lower right
4. Lateral teeth upper right & lower right


Install & use

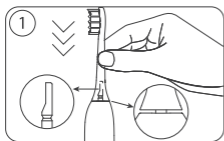
1. Put the head into the output axle carefully. Leave a space between the head and the main body, which is necessary for the vibration of the head. Try to start and test the toothbrush if working properly.

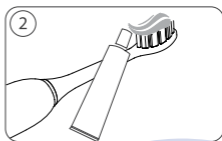
2. Put a little toothpaste on the head. Please don't turn on the toothbrush to prevent toothpaste from splattering.

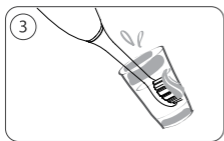3. Wet the head and toothpaste before brushing your teeth, which will allow the toothpaste to be more fully decomposed when brushing your teeth

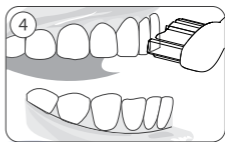4. Brush softly. Don't press too much, otherwise it will hurt gum and teeth surface. Move toothbrush head slowly, back and forth.

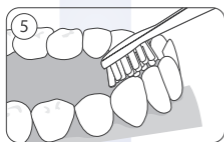5. Don't forcefully press teeth when brushing.

6. Rinse your mouth with water after brushing.

Charging instruction

■ Please start your toothbrush first when initial use, if you find the toothbrush is in the low power or no power situation, please make sure your toothbrush fully charged before use.

■ Without operation longer than 6 months, please charge your toothbrush fully according to the charging instruction (shown in the right pictures)


■ While using your toothbrush, if charging lamp keeps flashing and your toothbrush don't shakes.. Your toothbrush is in shortage of power. Please recharge in time.

■ Please put toothbrush body on the charger and confirm that the charger has safely powered on.

■ It indicates under the charging when the charge lamp keep red after putting the main body in charging base. And it indicates fully charged when the indicator light keep green flash.

■ With a long time non-recharging, the charging lamp may light up after 6-10 minutes when recharging. This is a normal situation.

- Charged as instruction will guarantee the service lifetime of the battery.
- The charging icon.

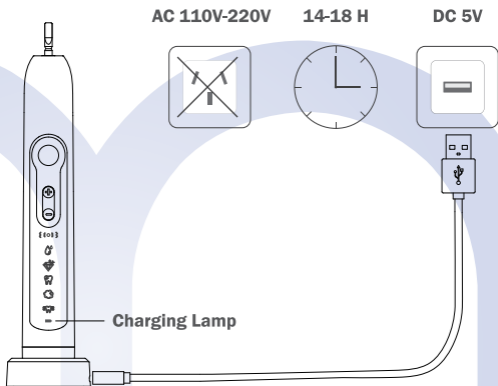Clean & Maintenance

Regular cleaning will maintain the best condition of the product and lengthen its service lifetime.

TO CLEAN THE BRUSH HEAD

Please clean the toothbrush head after every use, let it dry naturally.

Do not press and pull hard the bristle to avoid deformation and falling off.

TO CLEAN THE MAIN BODY

Clean your toothbrush with soft cloth and let it air drying. Do not immerse it in water for long time. Do not clean with hot water or acid. Cleaning solution: Suggest to use medical alcohol cleaning.

TO CLEAN CHARGER

Clean your toothbrush with soft cloth and let it air drying. Do not immerse it in water for long time. USB charging port do not immerse any liquid and do not clean with hot water or acid. Suggest to use medical alcohol cleaning.


TO REPLACE THE BRUSH HEAD

We suggest you change the head of your toothbrush every three months.

Environmental protection

■ •FL-A21 sonic toothbrush with a built-in rechargeable battery contains substances, which may pollute the environment. Please don't put this dumped toothbrush together with life garbage, need take out the battery and put it into recycling barrel. Our product should be delivered to corresponding recycling center, which is conducive to environmental protection. If you have any difficulties about taking out battery, you can send the product to our maintenance center.

■ If you intend to disassemble your product, please comply with the following procedures, and send all finished parts and battery to recycling center designated by the authority for environmental protection.

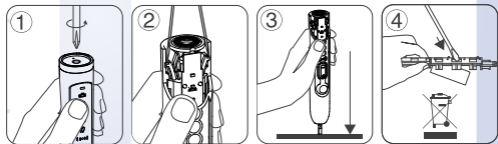❗ Discard rechargeable battery, this action is not reversible.

- According to the following steps, the operation shall comply with the basic safety precautions. Be sure to protect your eyes, hands and workbench surface.
- Remove toothbrush head from the main body, turning on the power until the battery is completely exhausted.
- Put the suitable crosshead(standard) screwdriver into the slot located at the base of the toothbrush holder, release the bottom cover by counterclockwise rotation of a screwdriver.
- Put the handle of your toothbrush on a hard surface with buttons upward, carefully open the enclosure from bottom edge with slot type screwdriver(standard), until the enclosure is 2.5 cm away the bed.
- The slotted screwdriver(standard) will be used to pry out the elastic buckle, meanwhile, pushing down the toothbrush shell, and taking out the inner component of the toothbrush.
- The slotted screwdriver(standard) will be inserted into plastic box to pry out the battery.

Fault diagnosis

■ TOOTHBRUSH WEEK VIBRATION

1. The toothbrush is forcefully pressed on the teeth. Press the brush head lightly on the teeth.
2. There is a foreign body between the head and toothbrush. Clean up the foreign body of head and toothbrush.
3. Did not adjust the Intensity to maximum strength. Intensity adjustment according to page 6.

■ TOOTHBRUSH CAN ONLY BE USED FOR A SHORT TIME EVEN IF IT IS FULLY CHARGED

1. The main body isn't fully charged in right way. Please charge the battery again, the main body is placed vertically on the top of the charger to make it contact with the charger. See page 11 for charging method.
2. The toothbrush is Insufficient charge time. Charging time have to be 12 hours approx.
3. The toothbrush is power off. Charging in time.

■ TOOTHBRUSH DOESN'T WORK

1. Toothbrush is power off. Please charge the battery again ,the main body is placed vertically on the top of the charger to make it contact with the charger. See page 11 for charging method.

2. Toothbrush with faults. Maintenance only at customer service center.

3. The toothbrush is reached its service lifetime. direction when there is a pause. (See page 8 for brushing teeth details). If toothbrush has been used for 3 years or more years, the declaration has reached the end of the period. Please timely return the processing according to the regulations

* If you do not understand, please consult the after-sales service center of manufacturer.

Technical parameters

Product name: ProBrush 9.0 (Ultrasonic)

Main Body

Rated voltage: 3.7V DC

Rated power: 1.7W

Charger

Rated voltage: 5.0V DC

Rated power: 2.3W

⚠ Please carefully read the parameters in the specification and use the product correctly.

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