Portable Oral irrigator

ProWater Clean 7.0

medica+

JAPAN TECHNOLOGY

Instruction Manual

Irrigator Station
ProWater Clean 7.0

Content

UL Standart	01 - 03
Product Description	04
Product Accessories	05
Operating Procedures	06 - 07
Usage Instructions	08 - 09
Recommended techniques	10 - 11
Disassembly Example	12
Battery Maintenance And Handling	13
Commons Problems	14 - 15
After-sales Services	16 - 17

UL Standart

IMPORTANT SAFEGUARDS

When using electrical products, especially when children are present, the basic safety precautions below should always be followed.

READ ALL INSTRUCTIONS BEFORE USING.

DANGERS:

To reduce the risk of electrocution:

- Do not handle charger with wet hands.
- Do not immerse in water or other liquid.
- Do not use while bathing.
- Do not place or store product where it can fall or be pulled into a tub, shower stall or sink.
- Do not reach for a product that has fallen into water. Unplug immediately.
- Check the charger cord for damage before the first use and during the life of the product.

Warning


To reduce the risk of burns, electric shock, fire, serious injury or damage to issue:

- Do not plug this device into a voltage system that is different from the voltage system specified on the device or charger.
- Do not use this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped into any liquid. Contact your local distributor.
- Do not direct water under the tongue, into the ear, nose or other delicate areas. This product is capable of producing pressures that may cause serious damage in these areas. See operating instructions for correct usage.
- Use this product only as indicated in these instructions as recommended by your dental professional.
- It's normal phenomenon that caused gum bleeding when you use the water flosser in the first time, which depended on different gum bearing ability.
- Fill reservoir with water or other dental professional recommended solutions only.


- Do not drop or insert any foreign object into any opening on the product.
- Keep charging unit and cord away from heated surfaces.
- Do not operate where oxygen or aerosol sprays are being used.
- Do not use mouthwashes/rinses that contain iodine, bleach, or tea tree oil.
- Do not while wearing any oral jewelry. Remove prior to use.
- Do not use if you have an open wound on your tongue or in your mouth.
- If your physician advised you to receive antibiotic premedication before dental procedures, you should consult your dentist before using this instrument or any other oral hygiene aid.
- For household use only. Do not use outdoors.
- Closely instruct and supervise children and individuals with special needs in the proper use of this product.
- Children should be supervised to ensure that they do not play with the appliance.

SAVE THESE INSTRUCTIONS.


Product Description

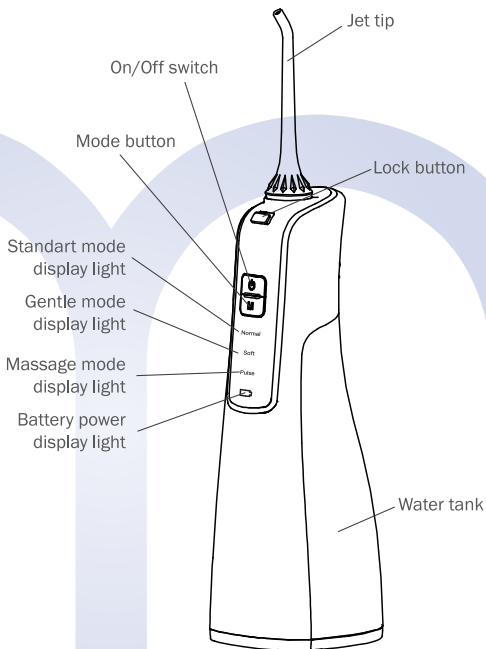Product Accessories

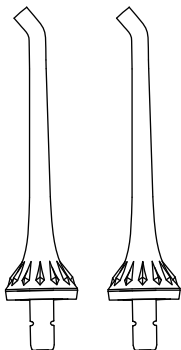Standart jet tips x 2 pcs.

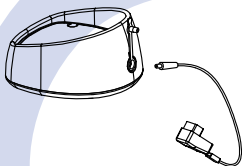Power recharger x 1 pcs.

The jet tip and other changeable components can be purchased at retail stores

Operator panel instruction


■ Preparation Work.

Extremely important! Prior to using for the first time, please fully charge the device, and insert the plug into the power socket, as well as connect to a power supply. During power charging the product cannot be used.

■ Installation And Dissambly of Jet Tip.

Insert the jet tip into the water outlet , until you hear a clicking sound.

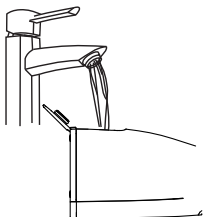As shown in Diagram 3, open the water tank lid and fill in water tank;

Note: please ensure that when filling in the water tank, the device is switched off.

If changing or removing the jet tip, please press on the disassembly button on the jet tip handle before removing the jet tip.


Usage Instructions

- This product has a memory function, and after switching it back on the water pressure mode will be the same as when the device was switched off.
- The standard, gentle and massage modes may be switched whether the device is on or off.
- The modes can be changed when the device is switched off, and you can select the mode you like before switching the device back on.
- Irrespective of the mode, after operating for 120 seconds the product will automatically switch off.
- When the switch button or mode button is used for with different colors, remaining battery power will be displayed. The red light show the battery power to be 35% at least. The blue light show the battery power to be 85% at least. The green light show the battery power to be full.
- During recharge, the three display lights will flash. When all three display lights are constantly on, this indicates that the power is recharged.

■ During usage, If all three display lights are flashing, this indicates that the product's power is insufficient, and recharging is required.


■ During recharge this product cannot be used.

■ Rotate the jet tip dial in order to adjust the direction of the water flow.

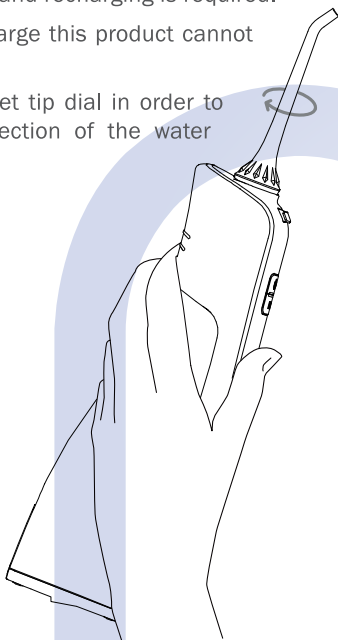
Recommended Techniques

Holding the device at downward slant, position the jet tip in the mouth, and direct the jet tip at teeth before switching on the device.

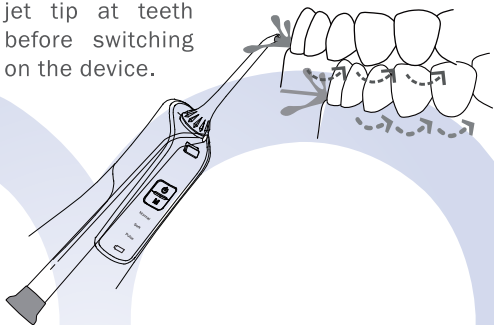
When the water flow and the gums are at a 90 degree angle, slightly close your mouth in order to avoid the water spraying out, and achieving the best results. Please first clench the teeth (rear teeth) and begin to apply to the incisors. Slide the water flow along the gums, as well as pause briefly between teeth, until the area around the gums and the spaces between the teeth are fully cleaned.

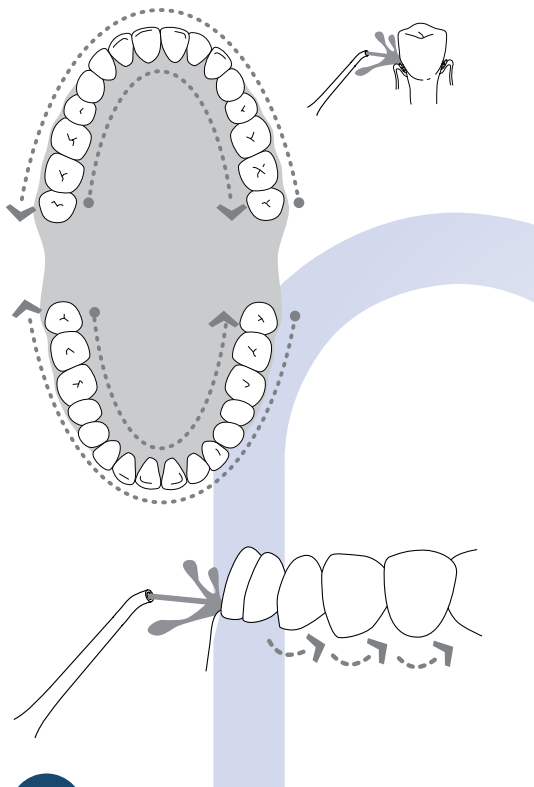
Disassembly Example

Holding the device at downward slant, position the jet tip in the mouth, and direct the jet tip at teeth before switching on the device.

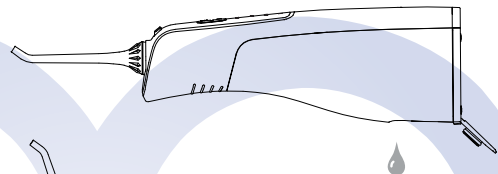2

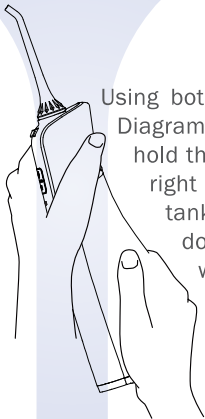Using both hands as indicated in Diagram 2 -use your left hand to hold the product body, and your right hand to hold the water tank. At the same time pull downwards to remove the water tank.

Battery Maintenance and Handling

Battery Maintenance

- Prior to first use the device it must be fully charged for approximately ten hours;
- If to be stored for a long-period of time, please fully charge the device in advance.
- In order to maximize the usage life of the battery, please promptly recharge when power is insufficient.

Battery Handling

- The device includes a lithium battery which cannot be changed. Towards the end of the battery life, please dispose in a local recycling centre.

Common problems and solutions

Problems	Reason	Method of handling
Leakage of water between jet tip and screw part	The seal ring on the fixed socket for the jet tip is damaged, and the jet tip cannot be installed properly	Contact the retailer for repair, and install the jet tip properly
Water tank leaks water	Seal ring is loose	Re-install correctly
Inadequate pressure	Low voltage usage	Use after recharging
Device can't be switched on	Low voltage usage	Use after recharging
Product is not discharging water	Filtration net is clogged	Clean the filtration net
Water tank sides are leaking water	Air entry hole is producing water	Standart phenomenon – high air pressure principle

If the above methods cannot resolve these problems, please consult with the retailer or proceed to the designated after-sales repair site for repairs.

Gentle reminder:

1. At the time of initial usage, please select “gentle mode,” in order to avoid the gums bleeding due to lack of adaption to usage. After consistent usage for two weeks there should be improvements. If bleeding is excessive or bleeding continues after two weeks, chronic gum inflammation could be an issue, and it is recommended that professional dental treatment be sought.
2. It is not recommended that the product be used by children, or people with dental implants, oral ulcers or open oral sores.

After-sales Services

Thank You For Purchasing The Medica+ Portable Oral Irrigator — Please Take Note Of The Following Issues:

1. Please retain the effective sales invoice issued by the retailer as a maintenance document.

2. Two years product malfunction warranty for non-human malfunctions is provided for this product in accordance with the relevant state provisions. If the product was purchased via Amazon, the documentation is not required, and only the order form number needs to be provided. In order for us to provide corresponding after-sales. Free maintenance (based on the invoice date for purchase of the device);

3. Where the following circumstances arise during the maintenance period,

costs and services fees shall be charged based on circumstances;

- The customer causes damage due to inappropriate usage or maintenance;
- Repairers who are not covered by third-party repair cause damage during repairs;
- Lack of a relevant sales document or more than two years have elapsed since the date of manufacture;
- Damage caused during usage for commercial purposes;
- Damage caused by force majeure events (such as storms, earthquake or fire).

4. Please carefully fill in the following form and deliver the device to the company's designated after-sale repairer for repair work.

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