## Photoepilator

**Hair Cleaner 8.0**

medica+

JAPAN TECHNOLOGY

medica+

JAPAN TECHNOLOGY

# **Instruction Manual**

## **Photoepilator**

**Hair Cleaner 8.0**

# Content

---

Introduction .....	01
Epilator Overview .....	03
How to use .....	05
Reference table .....	07
Safety information .....	09
What to Expect with the Photoepilator Hair Cleaner 8.0 Instrument IPL System? .....	17
Skin testing .....	19

## Introduction

Congratulations on the purchase of your new photoepilator Hair Cleaner 8.0 Instrument Intense Pulsed Light (IPL). Hair Removal system You're on your way for being liberated from the daily chore of hair removal.

the Photoepilator Hair Cleaner 8.0 Instrument system is a revolutionary light-based device designed for the removal of unwanted body hair at home.

It uses the same Intense Pulsed Light (IPL) technology that professional derma-otologists and salons use, but is optimized for use for/of the privacy and comfort of your own home.

Unlike most at-home hair removal methods, the Photoepilator Hair Cleaner 8.0 Instrument system doesn't cut or pull out hair-it actually works below the skin's surface to disable active hairs and helps prevent new ones from growing.

This is a gradual process. You'll start to see the results of two weeks after your initial treatment and results will improve gradually on time.

How does the Photoepilator Hair Cleaner 8.0 Instrument System Work?

The process of light-based hair removal is well-known and has been clinically prov-en around the world for over 15 years. the Photoepilator Hair

Cleaner 8.0 Instrument system uses the same proven technology to emit an intense pulse of light (IPL) to heat up and disable active hair follicles.

The light energy targets dark coloring in body hair, which is why it works best on lighter skin tones and naturally dark hair.

To determine if your skin and hair color are ok to use with the Photoepilator Hair Cleaner 8.0 Instrument system, refer to information provided. Every hair in our body goes through a continuous growth cycle, from active to resting phases. The time it takes to complete a full hair growth cycle varies, but is typically 18-24 months. Light-based treatments only affect hairs in the active phase, which can account for to 90% of hair follicles.

That's why we recommend using the Photoepilator Hair Cleaner 8.0 Instrument system.

Three times in the first month to ensure you hit every active hair.

After that, use the Photoepilator Hair Cleaner 8.0 Instrument system as needed for any missed hairs or hairs that have recently become active.

■ The first 3 sessions with the Photoepilator Hair Cleaner 8.0 Instrument IPL system should be 2 weeks apart.



■ After that, use the Photoepilator Hair Cleaner 8.0 Instrument system when you start to see hair re-growth in the treated area.

■ Do NOT treat the same body part more often than every 2 week. After a hair follicle is treated, it can take one to two weeks for the hair to fall out. During this time, it may look like the hair is still growing, but really the hair is being "pushed out". You may or may not see these stubs pushed out from the skin.

If you start treating a different body part, start the process of for that body part.

## **Photoepilator Hair Cleaner 8.0 instrument system diagram reference**

1. Hair Cleaner 8.0 Instrument,
2. Screen,
3. Skin contacting area,
4. Flash button,
5. Power switch
6. Socket for power,
7. Flash display (digital) to the remained service
8. Gear indicator,
9. Auto continuous flash mode,
10. Fan
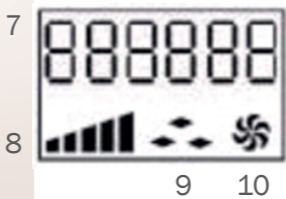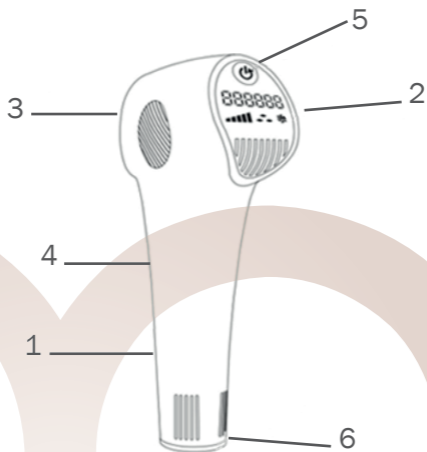
# How to use your photoepilator Hair Cleaner 8.0 instrument system

## HOW TO USE

1. Please remove the hair on the skin surface firstly. If the hair does not exceed 2 mms, the instrument can be operated directly on shaving.
2. Please clean the skin and make sure it is hair-free, dry and free from oily substances. It is strictly prohibited to apply condensation, aloe vera gel and other products.

## STEPS

1. Please clean the skin surface and filter with a soft cloth before use.
2. Please connect the power adapter and intensity indicator. The flashing " red light" indicates that the circuit is normal.
3. Please press the power switch until the screen lights up.
- 4 Put the device at a 90° angle on the skin and make sure the flash window is full contact with your skin.
5. If the screen changes from "flashing " onto Constantly bright", It indicates the depilator is electrified.
6. After the machine is pressed on the skin, please gently press the flash button until it emits a flash.

7. After using the instrument, please replace the depilated part.

8. Please shine as close as possible to the site of previous treatment to ensure that all sites are adequately treated.

One hair removal treatment is performed 3 times. After flashing the position on first hair removal, then do the second and third times, do not continue to irradiate 3 times for/of the same place.

### **MODES**

Two treatment modes including [press and flash] and [slide and flash] are convenient for you to use on different body parts:

1. The press and flash mode is suitable for small or curved areas, such as knees and armpits. You just need to press the flash button, then the flash will appear.

2. The slide and flash mode is convenient for large areas such as legs. Please hold down the flash and slide the instrument along the skin at one side to emit multiple flashes continuously.

### **WHAT IS THE INTENDED USE FOR THE PHOTO EPILATOR HAIR CLEANER 8.0 INSTRUMENT IPL SYSTEM?**

the Photoepilator Hair Cleaner 8.0 Instrument system is an over-the-counter device intended for

the removal of unwanted hair.


**Skin Tone** It is intended for use only on light to medium skin colors.

To determine if your Skin Tone is suitable to use with Hair Cleaner 8.0 Instrument places skin in each intended treatment area next to the Skin Tone Chart.

If your Skin Tone matches panels 1-4, You can use the Photoepilator Hair Cleaner 8.0 Instrument system. If your Skin Tone matches 5-6, DO NOT use the system.

Always test each different area you wish to treat with the Skin Tone Chart and Skin Tone Sensor before treating. Some areas are darker than others.

## Theraupctic cycle reference table
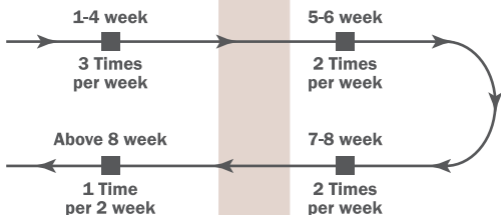
## HAIR COLOR

It is intended for use for naturally black or brown hair. Do not use Hair Cleaner 8.0 Instrument on white or gray hair because it will not work.

IPL is less effective against lighter colors, such as blond or red hair.

**Body Areas** It is intended for use by women and men on body hair below the neck, including legs, underarms, bikini line, chest, stomach, and arms.

It has NOT been cleared in the U.S. for use for the face, head, ears, back or neck.

## Safety information

Please read and consider the information on the following section before use.

DO NOT use photoepilator Hair Cleaner 8.0 Instrument on Dark Skin.

Darker skin may absorb too much light energy which could cause skin swelling, temporary blisters, or could change the color of your skin.

if skin tone in the intended treatment area is suitable for use for the Photoepilator Hair Cleaner 8.0 Instrument system.

You should test each different area of skin you want to treat.

Some areas are darker than others.

DO NOT use on tattoos, dark brown or black spots (such as dark freckles, birthmarks, moles or warts). Tattoos or dark spots may absorb too much light energy which could cause skin swelling, temporary blisters or could change the color of your skin .

Refer to the chart to determine if tattoos or dark spots in the intended treatment area are suitable for use with Hair Cleaner 8.0 Instrument.

You should test each different area of skin you

want to treat. Some areas are darker than others. DO NOT use on the face, head, ears, neck, nipples, genitalia, or around the anus. These areas may have darker skin and using the Photoepilator Hair Cleaner 8.0 Instrument there could cause skin swelling, temporary blisters, or could change the color of your skin.

DO NOT use on or around the eyes, eyebrows or eyelashes.

Doing so can cause serious and permanent eye injury.

## **WARNINGS**

the Photoepilator Hair Cleaner 8.0 Instrument system is a powerful device .

As such it should be used with special attention to safety

Please read all warnings and safety precautions before use and strictly follow them when using the Photoepilator Hair Cleaner 8.0 Instrument system .

DO NOT flash the Photoepilator Hair Cleaner 8.0 Instrument on the identical treatment spot multiple times of the same session, as this could cause skin swelling, temporary blisters or could change the color of your skin .

DO NOT use on damaged skin, including open wounds, cold sores, sunburns, or areas where

you have had a skin peel or other skin resurfacing procedures within the last 6-8 weeks (check with your doctor).

DO NOT use on irritated skin, including rashes or swollen skin.

These conditions can make the skin more sensitive to the light treatment.

Wait for the affected area to heal before using .

DO NOT use if your skin is sensitive to light.

In order to determine if you can use, perform a test of a small patch of skin and wait 24 hours to determine if there are any adverse reactions before using more broadly.

DO NOT use if you are or may be pregnant, or are breast-feeding.

the Photoepilator Hair Cleaner 8.0 Instrument IPL system has not been tested for pregnant women.

Hormonal changes could increase skin sensitivity and the risk of skin injury.

DO NOT use for any purpose other than hair removal.

DO NOT flash the Photoepilator Hair Cleaner 8.0 Instrument against any material other than skin.

DO NOT flash into an environment where a spark could cause an explosion, such as near a gas

leak. DO NOT use on any area where you may someday want your hair to grow back.

The results may be permanent and irreversible.

DO NOT use the system if it appears damaged, and DO NOT attempt to open or repair the Photoepilator Hair Cleaner 8.0 Instrument system .

These actions increase the risk of skin or eye injury to the flash cartridge. Contact Customer Service for assistance.

DO NOT allow the hand piece or base to become damp or wet.

DO NOT use near bathtubs, showers, basins, or other vessels containing water.

DO NOT operate in a damaged plug or cord .

Keep the cord away from heated surfaces (such as stoves, fireplaces, radiators or furnaces).

DO NOT use if any part of the Photoepilator Hair Cleaner 8.0 Instrument system is cracked or broken. Contact Customer Service for assistance.

DO NOT open the unit (except to replace the CARTRIDGE ) or attempt to repair your system as there are no user serviceable parts of the Photoepilator Hair Cleaner 8.0 Instrument system.

Trying to open it may also damage the system and will void your warranty. Please contact

Customer Service if you have a broken or damaged system in need of repair.

Failure to follow these precautions can cause electrocution or electric shock .

DO NOT cover any part of the Photoepilator Hair Cleaner 8.0 Instrument BASE UNIT or HAND PIECE ventilation slots into,as they help to keep the system cool during operation.

Covering the ventilation slots could result in over-heating.

KEEP OUT OF REACH OF CHILDREN the Photoepilator Hair Cleaner 8.0 Instrument IPL system is a powerful device and should not be used by anyone under the age of 18.

DO NOT leave the system unattended around children,as the cords represent a choking hazard. Children should be supervised to ensure that they do not play with the appliance.

### **BENEFITS AND RISKS.**

A medically-supervised clinical study showed that the Photoepilator Hair Cleaner 8.0 Instrument IPL system was effective against the removal of unwanted hair when used as directed.

60 participants entered the clinical trial with the Photoepilator Hair Cleaner 8.0 Instrument and treated 3 times of a 4 week period.

Treatments were performed on underarms, legs, bikini area, stomach, chest and arms. Hair counts were performed on site before treatment and again 3 and 6 months after the final treatment. Any side effects observed or reported were documented by medical personnel. This study showed the benefits and risks of hair removal of the Photoepilator Hair Cleaner 8.0 Instrument system . The following sites are contraindicated and were therefore excluded from treatment: face, Head, ears, neck, nipple, genitalia, anus, eyes, eyebrows or eye lashes.

**Benefits** When used as directed, Hair Cleaner 8.0 Instrument IPL system helps remove unwanted hair without the need for daily shaving, wax- ing, or plucking . In a clinical study, after just three treatments, participants show An average of 37% less hair 3-months after the last treatment and 48% less hair 6-months after last treatment. Actual results vary from person to person, depend on the body area being treated, and whether instructions were followed correctly.

**Risks** It is important to use the Photoepilator Hair Cleaner 8.0 Instrument IPL system properly. Make sure that you read and follow all the warnings and that you don't use the device if it

isn't right for you . If you don't understand any of the warn- ings or have any questions, please contact Customer Service .

If you use the Photoepilator Hair Cleaner 8.0 Instrument system as directed,the risks are low. In a clinical study,most users reported mild or no discomfort at all from treatmines,with sensation decreasing over multiple treatment sessions.

This is normal. You may feel warmth,burn- ing,tingling,or itching; some users described feeling the same level of discomfort as when a rubber band is lightly snapped against the skin. If the pain is too intense or persists after a treatment, stop using and contact Customer Service. Your skin may become red within 24 hours after treatment.

This redness generally clears up within a few hours. See your doctor if the redness does not go away within 2-3 days

Temporary blisters .These are unlikely to people with a light or medium skin tone .

If you do develop a blister,stop using Hair Cleaner 8.0 Instrument, clean the blister with soap or antiseptic and cover the blister with a bandage and change daily to help prevent the blister from becoming infected. Contact Customer Service if you need more information.

In rare cases, some users may see an increase in the number of hairs in the treated area after treatment. Since hair grows in cycles, this increase is generally temporary .

If at any time after you begin using Hair Cleaner 8.0 Instrument you believe that the number of hairs regrowing in the treatment area is greater than the number that existed before you started treatment, stop using Hair Cleaner 8.0 Instrument and contact Customer Service.

If this increase persists with several treatments, stop using and see your doctor.

The following symptoms affected no participants in the clinical study, but have been reported on / to other hair removal lasers or IPL devices .

- Skin swelling, especially with the hair follicle. If the swelling doesn't go away within 2-3 days, stop using and contact Customer Service .

- Skin discoloration, If you notice that your Skin Tone in the treatment area has changed, stop using and contact Customer Service.

- Eye injury. Flashing any laser or IPL device near or around the eyes could cause serious and permanent eye injury.

This is very unlikely so long as the device is used below the neck as indicated .

## WHAT TO EXPECT WITH THE PHOTOEPILATOR HAIR CLEANER 8.0 INSTRUMENT IPL SYSTEM?

During Treatment, the Photoepilator Hair Cleaner 8.0 Instrument IPL system has been designed for ease of use, and hair removal sessions should go by quickly. During a session it is normal to experience and feel:

- **A Fan Noise:** the fan starts up once the Skin Tone Sensor is unlocked. The noise is similar to, but quieter than a hairdryer. This is normal and is key to keeping the unit and flash lamp cool.

- **A Flash of Light:** the Skin Contact Sensors ensures that most of the light is directed to the skin. The remaining visible light will not harm the eyes when applied to non-facial sites, and special eye protection is not needed as long as you follow directions in this Use and Care Guide.

- **A Sound with Each Flash :** When a flash of light is activated, it is normal to hear a subtle "pop" sounds as the light flashes.

- **A Sensation of Warmth and Tingling :** Actual sensations felt by users will vary by person and body part being treated. Many people describe each flash like a rubber band being lightly snapped against the skin. If the sensation is too intense or painful, stop using and contact Customer Service.

■ **Some Mild Redness** : During and just after your IPL session, it is not uncommon to see some very mild redness, mostly noticeable around the hairs themselves. This is normal and typically clears up within a few hours. However, contact Customer Service if the redness does not go away within 2-3 days.

### **THERAPEUTIC CYCLE REFERENCE TABLE**

One hair removal treatment is performed for 3 times .

After flash the position on first hair removal, then do the second and third times,do not continue to irradiate 3 times for the same place .

The higher the gear position , the better the effect.

In the early stage, it is recommended to using 1-2 gears first After the skin is adapted , you can choose to use gear 5 .

The interval between treatments can be determined according to your personal hair growth rate .

It will be different from one another to the locations of your body .

Skin Test Before using for the first time or after using it after sunburn, perform a skin test of each treatment area.

## **SKIN TESTING IS A NECESSARY STEP TO CHECK THE SKIN'S RESPONSE TO THE TO THE COURSE OF TREATMENT AND TO DETERMINE THE APPROPRIATE LIGHT INTENSITY POSITION FOR EACH BODY PART.**

1. Select the area close to the course of treatment.
2. Turn on the instrument and confirm that gear 1 is selected.
3. Press the instrument against the skin and press the flash button to give a flash .
4. The instrument slides into the skin to the next hair removal position.
5. Adjust the gear up one level and send a flash, and then slide the instrument to the next position.
6. After completing the skin test, please wait 24 hours to check if there is any reaction to the skin. If the skin responds. choose the most suitable gear that does not cause skin reactions and use it in the subsequent hair removal sessions.

### **AFTER USE**

- 1 .The skin may have a slight redness or aching or warm feeling .  
These reactions are completely harmless and will soon disappear.
2. You can safely apply lotion,cream,body

lotion,moisturizer or cosmetics to the treatment area . If the skin is stinging or red and swollen after the end of the treatment;please wait until these conditions disappear and apply any product of your skin .

If the skin has a tingling sensation after applying the product,wash it off with water.

		SKIN TONE					
							
HAIR COLOUR		✓	✓	✓	✓	✓	×
		✓	✓	✓	✓	✓	×
		✓	✓	✓	✓	✓	×
		×	×	×	×	×	×
		×	×	×	×	×	×
		×	×	×	×	×	×

✓ = suitable                      × = not suitable / effective for use

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Treatments were performed on underarms, legs, bikini area, stomach, chest and arms. Hair counts were performed on site before treatment and again 3 and 6 months after the final treatment. Any side effects observed or reported were documented by medical personnel. This study showed the benefits and risks of hair removal of the Photoepilator Hair Cleaner 8.0 Instrument system. The following sites are contraindicated and were therefore excluded from treatment: face, Head, ears, neck, nipple, genitalia, anus, eyes, eyebrows or eye lashes. Benefits When used as directed, Hair Cleaner 8.0 Instrument IPL system helps remove unwanted hair without the need for daily shaving, waxing, or plucking. In a clinical study, after just three treatments, participants show an average of 37% less hair 3-months after the last treatment and 48% less hair 6-months after last

treatment. Actual results vary from person to person, depend on the body area being treated, and whether instructions were followed correctly.

**Risks** It is important to use the Photoepilator Hair Cleaner 8.0 Instrument IPL system properly. Make sure that you read and follow all the warnings and that you don't use the device if it isn't right for you . If you don't understand any of the warnings or have any questions, please contact Customer Service.

If you use the Photoepilator Hair Cleaner 8.0 Instrument system as directed, the risks are low. In a clinical study, most users reported mild or no discomfort at all from treatments, with sensation decreasing over multiple treatment sessions.

This is normal. You may feel warmth, burning, tingling, or itching; some users described feeling the same level of discomfort as when a rubber band is lightly snapped against the skin. If the pain is too intense or persists after a treatment, stop using and contact Customer Service. Your skin may become red within 24 hours after treatment.

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**Temporary blisters.** These are unlikely to people with a light or medium skin tone. If you do develop a blister, stop using Hair Cleaner 8.0 Instrument, clean the blister with soap or antiseptic and cover the blister with a

bandage and change daily to help prevent the blister from becoming infected. Contact Customer Service if you need more information. In rare cases, some users may see an increase in the number of hairs in the treated area after treatment. Since hair grows in cycles, this increase is generally temporary .

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