ENG / UA / RU

Kids Brush 2.0

medica+
JAPAN TECHNOLOGY

medica+

JAPAN TECHNOLOGY

ENG

Instruction manual


Kids Brush 2.0

Content


Use Steps	01 - 02
Parameters	03
Functions.....	04-05
Characteristics	06-07
Product Structure	07
Attentions	08
Warranty terms	09-10

Use Steps


First use

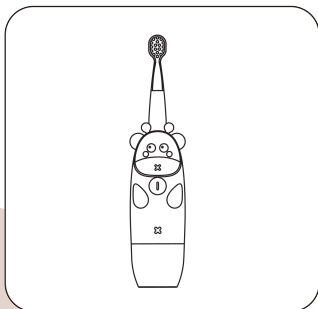
1. Take out the toothbrush.

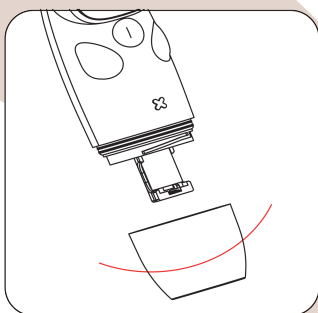2. Anticlockwise turn the bottom cover and open the bottom cover.

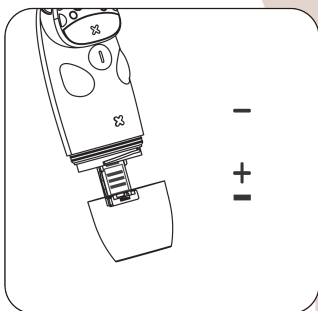
3. Take out the battery and the battery protective film, and replace the battery (Positive down).

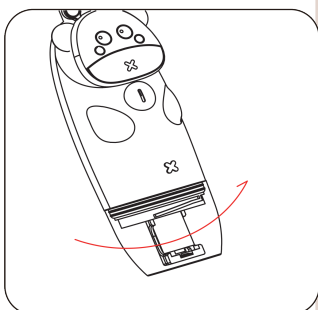
4. Cover the bottom cover clockwise screw it down.

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
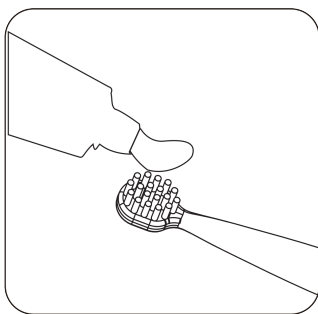
5. Squeeze toothpaste on the brush.

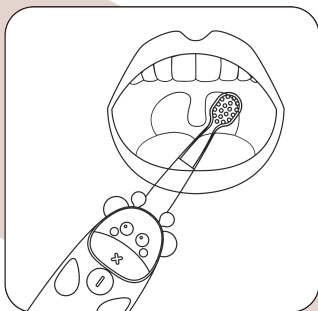6. Put the toothbrush into your mouth.

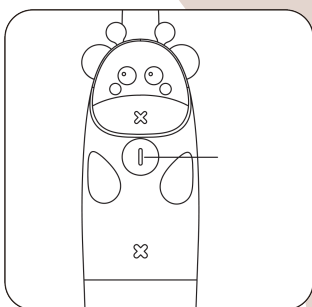
7. To select an appropriate gear, gently press the switch. It is 2 minutes for timing brush, and there is a reminder for changing the brushing area in 30 seconds. Press and hold for 3s to turn it off.

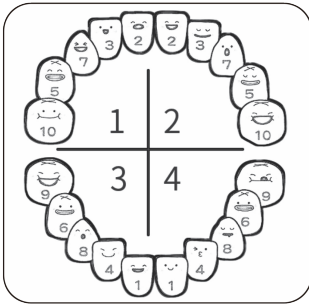8. Brush scientifically: each tooth area shall be brushed for 30s for cleaner teeth.

9. Rinse the mouth with water cup and clean toothbrush with water after brushing the teeth.

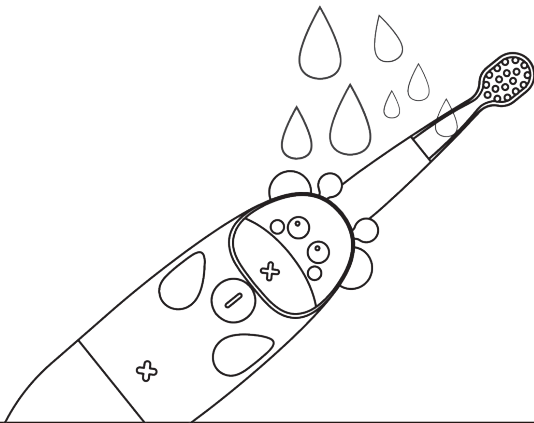Parameters

Main Engine

Product Size: 160.5 * 35 * 35mm

Rated Voltage: 1.5V

Rated Power: 0.51W

Rated Current: 0.34A

Functions

1. Four-gear Mode

First Gear Seven-color Mode: turn it on and display red, orange, yellow, green, blue, indigo and purple, the interval of color change is 0.5 without vibration, it lasts for 1 minute. This gear is recommended for children who are not suitable for the earthquake.

Second Gear Gently Brushing Mode: keep the light unchanged, it lasts for 2 minutes. This gear is recommended for children aged 2-5 years.

Third Gear Standard Brushing Mode: the light changes to the next color after the second gear, it lasts for 2 minutes. This gear is recommended for children aged 5-8 years.

Fourth Gear Strong Brushing Mode: the light changes to the next color after the third gear, and the total time is 2 minutes. This gear is recommended for children aged 8-12 years.

2. 2-minute Timing

For the second, third and fourth gear modes, when it is 30s, 60s and 90s, the lights flash once, twice and three times respectively, and the motor stops for one time, two times and three times respectively. The flashing and stopping time is 0.5s to remind changing of brushing area. 120s, the light flashes twice, and then it will be turned off.

3. Memory Function


Under any state of use, it will enter the memory mode after long press to turn off or shut down automatically. It will automatically enter the last use mode when you start the toothbrush next time.

Characteristics

DuPont Brush

It is made of American DuPont material, the end is rounding, which is strong and wear-resistant. It can protect the baby's gums while cleaning greatly.

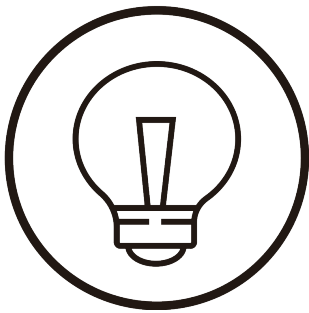Colorful Lights


It has cute animal model with colorful glare, children will love it for a good habit of teeth brushing.

Silicone Sleeve Design

The unique and creative silicone design is three-dimensional with more comfortable and smoother feeling, and kids will love the feeling of teeth brushing.

Product Structure

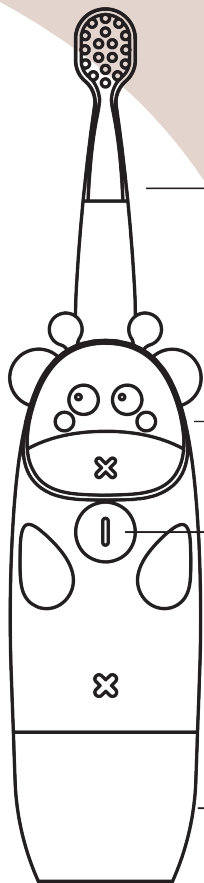Replaceable Toothbrush Head

Main Engine

On-off Key

Battery Bottom Cover

Attentions

1. The normal working temperature is 0-50°C.
2. It is not used for a long time, please remove the battery, and store the toothbrush in a box and place it in a cool place.
3. Please try your best to dry the toothbrush head after brushing teeth.
4. When the product is working, the output power is maximum. It is forbidden to let the product work continuously for a long time.
Generally, the use time should be less than 10 minutes continuously.
5. Never walk or run when putting the toothbrush in your mouth.
6. Make sure that a supervisor is accompanying with kids and never leave your sight when brushing teeth.

Warranty terms

Dear customers:

Thanks for your support for the brand of MEDICA+. Please keep this card properly and enjoy the warranty service-with the warranty card.

Since the date of purchase, the warranty is guaranteed for one year in case of failure under normal conditions!

The following conditions are not covered by the warranty.

1. The product is damaged due to improper personal use or lack of reasonable maintenance;
2. The product is damaged by dismantling and reforming without permission;
3. The parts which do not belong to our company are applied to the product;
4. Damage caused by force majeure.

For the paid maintenance service provided outside the warranty, charge the corresponding fee according to the regulations. When applying for warranty, please send the product together with the documents (for example, invoice) which can prove the date of purchase to the customer service department.

The parts or the whole machine replaced during the warranty period belong to MEDICA+.

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