
ENG / UA / PL


**Portable Oral
irrigator**

ProWater Clean 7.0

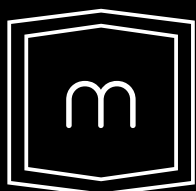medica+
JAPAN TECHNOLOGY

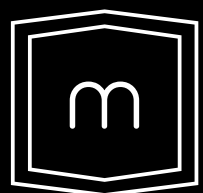
medica+

JAPAN TECHNOLOGY

ENG

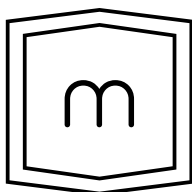medica+
JAPAN TECHNOLOGY

Instruction manual

ProWater Clean 7.0

Contents

Caution.....	01-03
Product Description.....	04
Product Components.....	05
Control Panel.....	06-07
Usage.....	08-09
Recommended Methods.....	10
Disassembly.....	11
Battery.....	12
Common Issues.....	13-14
Warranty.....	15
Technical Specifications.....	16

Dear Customer!

Thank you for choosing Medica+ products.

Please read this manual carefully before using the device.

Caution

When using the device, especially near children, it is important to follow basic safety precautions, including:

To reduce the risk of electric shock:


- Do not touch the device with wet hands.
- Do not immerse it in water or any other liquid.
- Do not use it while bathing.
- Do not place or store the device in areas where it may fall into a bathtub, shower, or sink.
- Check the charging cord for damage before first use and throughout the device's lifespan

To reduce the risk of burns, electric shock, fire, serious injury, or device damage:

1. Do not use the device if the cord is damaged, if it is not functioning properly, or if it has been immersed in liquid. Contact a service center.
2. Do not direct water under the tongue, into the ear, nose, or other sensitive areas. The device may generate pressure that could cause serious injury.
3. Use this device only as specified in the manual and as recommended by your dentist.
4. Bleeding gums during the first use of the irrigator is normal and depends on individual characteristics.
5. Fill the reservoir only with water or other dentist-approved solutions.
6. Do not drop or insert foreign objects into the device's openings.
7. Keep the charger and cord away from heated surfaces.
8. Do not use the device in areas where oxygen or aerosol sprays are being used.
9. Do not pour mouth rinses containing iodine, bleach, or tea tree oil into the reservoir.


10. Do not use the device while wearing oral jewelry. Remove it before use.
11. Do not use the device if you have an open wound on your tongue or in your mouth.
12. For household use only. Do not use outdoors.
13. Instruct and supervise children and individuals with special needs on the proper use of this device.
14. Children should be supervised to ensure they do not play with the device.

Product Description

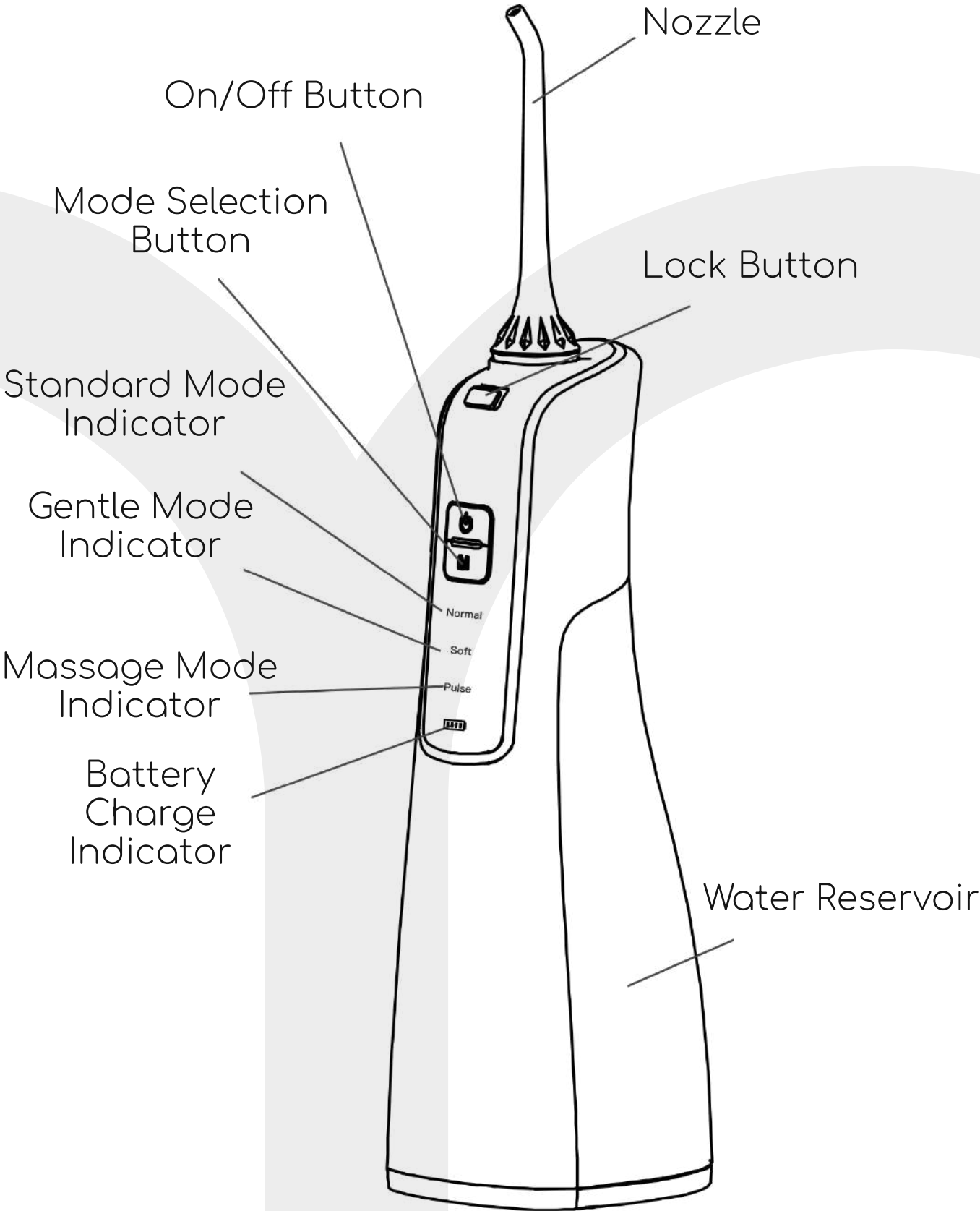
Product Components

- Irrigator – 1 piece
- Nozzles (0.8 mm) – 2 pieces
- Warranty Card – 1 piece
- User Manual – 1 piece

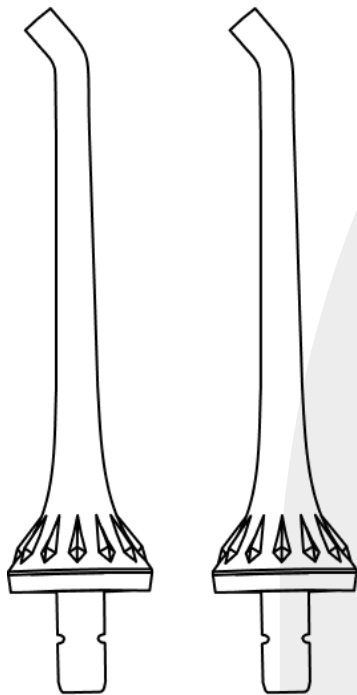
Standard Nozzles
x 2 pieces

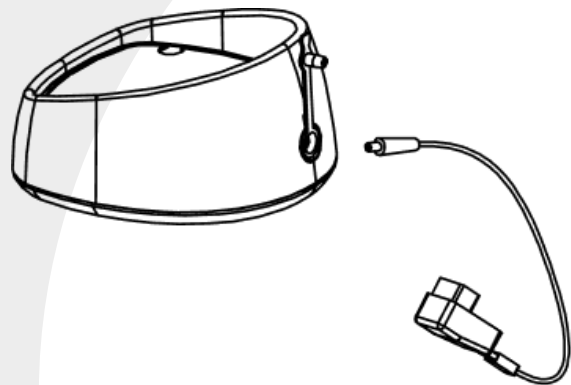
Charger x 1 piece

Control Panel

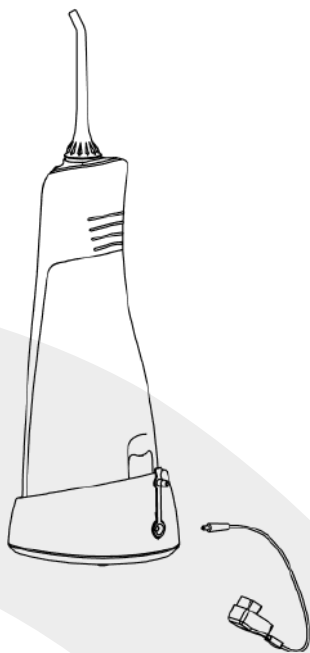
Preparation

Before first use, fully charge the device by plugging the adapter into a power outlet and connecting it to the device. The device cannot be used while charging.

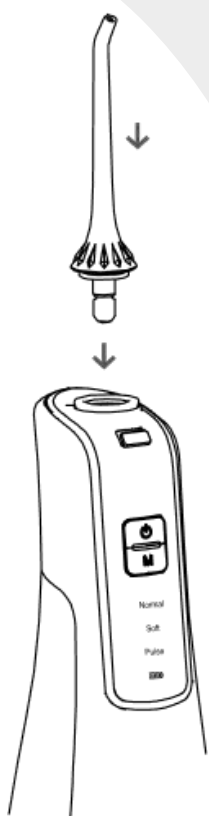
Installing and Removing Nozzles

Insert the nozzle into the water outlet until you hear a clicking sound. Open the water reservoir lid (as shown in the illustration) and fill it with water. The device must be turned off while filling the reservoir. To replace or remove the nozzle, first press the designated button, then gently pull it out.

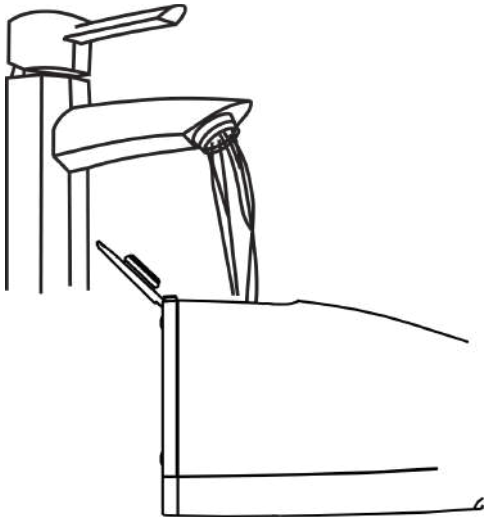
Open the water reservoir lid as shown in the illustration and fill it with water. Ensure the device is turned off while filling the reservoir.

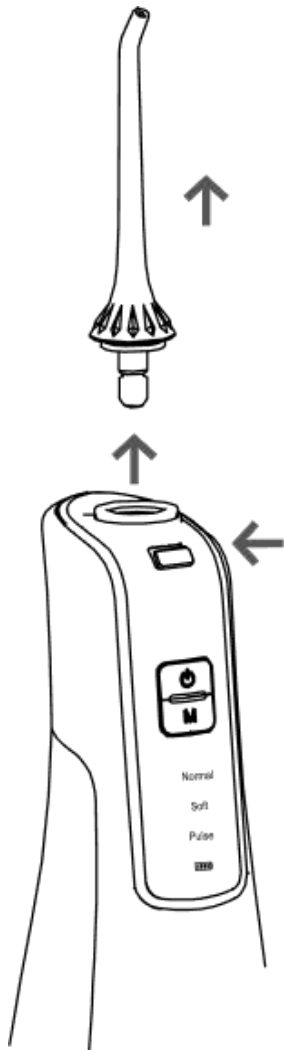To replace or remove the nozzle, press the designated button before pulling it out.

Usage

1. The device has a memory function, so when turned on again, the water pressure mode will remain the same as it was before turning off.
2. Normal, Soft, and Pulse modes can be switched whether the device is on or off.
3. Modes can be changed while the device is off, allowing you to select the desired mode before turning it on.
4. Regardless of the mode, the device will automatically turn off after 120 seconds of operation.
5. If the power button or mode selection button lights up in different colors, it indicates the battery level:
 - Red light: Battery charge is at least 35%.
 - Blue light: Battery charge is at least 85%.
 - Green light: Battery is fully charged.
6. If all three indicators flash simultaneously during use, the battery is low and needs to be recharged.
7. Rotate the nozzle to adjust the direction of the water flow.

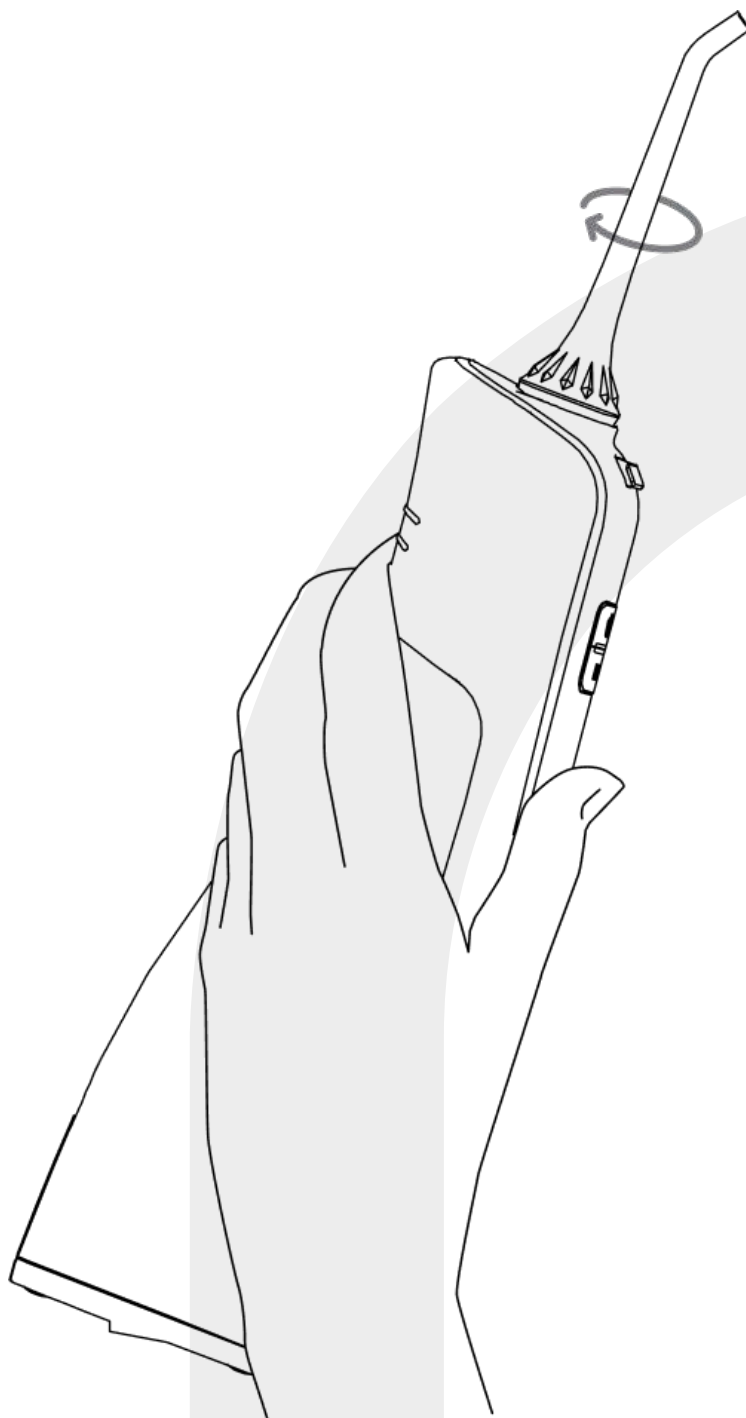
Recommended Methods

Hold the device at a downward angle, place the nozzle in your mouth, and aim it at your teeth before turning it on.

When the water stream and gums form a 90-degree angle, slightly close your mouth to prevent water splashing.


Begin by cleaning the back teeth, then move towards the front teeth (incisors).

Guide the water stream along the gumline, pausing briefly between teeth until the gum area and interdental spaces are thoroughly cleaned.

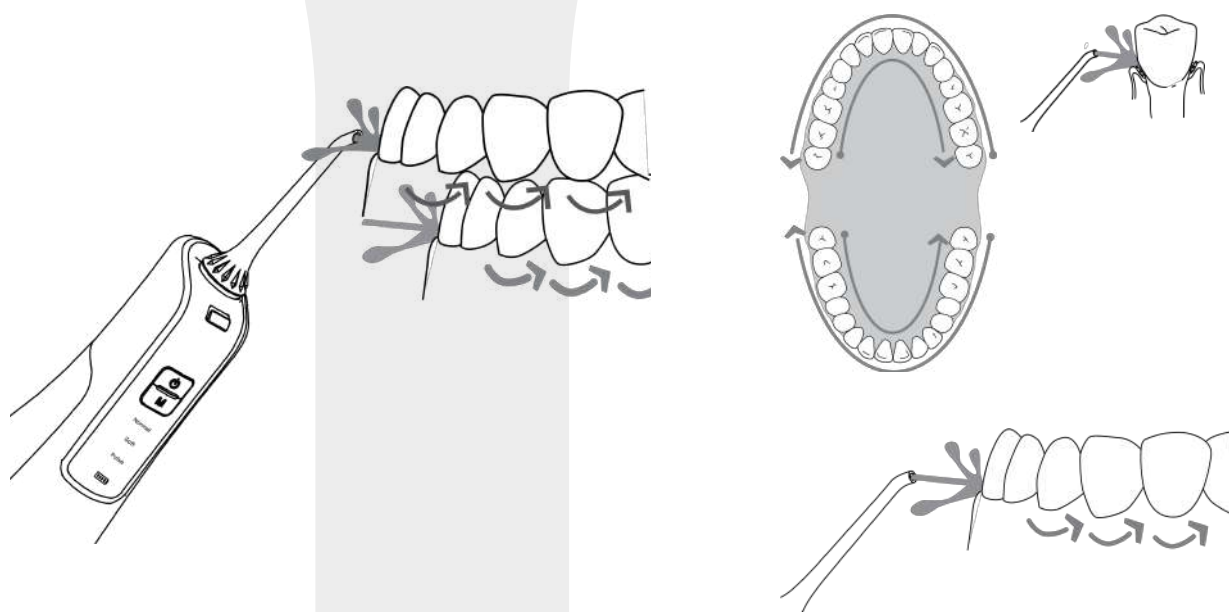Disassembly

To detach the water reservoir from the device, hold it with both hands as shown in the illustration. Hold the main body of the device with your left hand and the water reservoir with your right hand. Pull the reservoir downward to remove it.

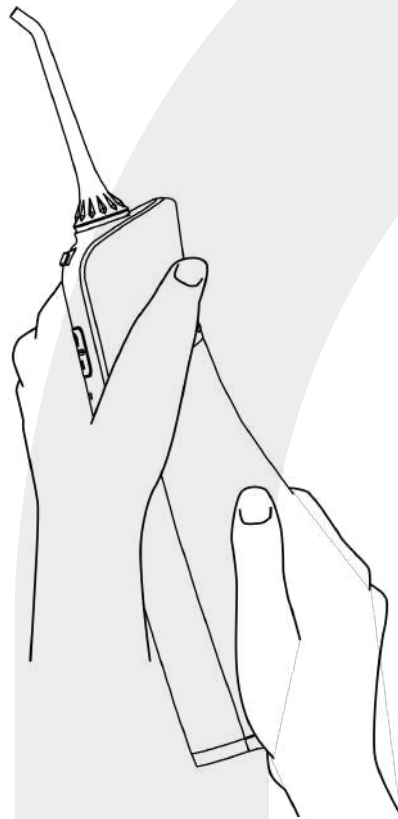Battery

Battery Maintenance

- Fully charge the device for ten hours before first use.
- If you do not plan to use the device for an extended period, fully charge it in advance.
- To maximize battery lifespan, charge it when the battery level is low.

Common Issues

Issue	Cause	Solution
Water leakage between the nozzle and the threaded part	The sealing ring on the fixed socket for the nozzle is damaged, and the nozzle cannot be properly installed.	Contact the service center.
The water reservoir is leaking.	The sealing ring is not tightly fitted.	Reinstall it correctly.
Insufficient pressure.	Low battery charge.	Use after recharging.
The device does not dispense water.	The filter screen is clogged.	Clean the filter screen.
Water is leaking from the sides of the reservoir.	Water is entering the air intake hole.	A standard phenomenon – the principle of high air pressure.

If the above methods do not resolve these issues, please call the service center. The phone number is listed in the warranty card.

Important Reminder:

1. When using the device for the first time, select the "Soft mode" to avoid gum bleeding due to insufficient adaptation. After consistent use for two weeks, improvement should occur. If bleeding is excessive or persists after two weeks, it may indicate chronic gum inflammation, and it is recommended to consult a specialist for dental treatment.
2. The product is not recommended for use by children, as well as individuals with dental implants or open sores in the oral cavity.

Warranty

Dear Customers!

Please keep this instruction manual properly to use the warranty certificate. The warranty is provided for one year from the date of purchase (or for 2 or 3 years if an extended warranty is purchased) in case of malfunction under normal operating conditions.

The warranty does not cover the following cases:

1. If the product is damaged due to improper use or lack of proper maintenance.
2. If the device is damaged due to disassembly or modifications.
3. If parts not belonging to our company have been used with the device.
4. Damage caused by force majeure circumstances.

According to the terms, maintenance services provided outside the warranty period are subject to additional charges.

Parts replaced during the warranty period remain the property of MEDICA+

Technical Specifications

Water tank capacity: 170mL

Input voltage: 100-240VAC 50/60Hz;

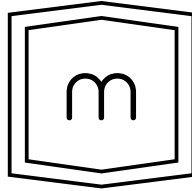
Battery: lithium battery with a capacity of 1400mAh

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