# Ultra face scrubber

VibroSkin 5.0

medica+

JAPAN TECHNOLOGY

medica+

JAPAN TECHNOLOGY

# **Instruction manual**

**VibroSkin 5.0**

## Contents

Description .....	01-03
Important Safety Measures .....	04
Indications .....	05
Contraindications .....	05-06
Usage Recommendations .....	06
Mode Descriptions .....	07-10
Technical Specifications .....	11
Warranty Terms .....	11-12
Package Contents .....	12

# Description

Dear Customer!

Thank you for choosing Medica+ products. Please read this manual carefully before using the scrubber.

Ultrasonic vibrations are high-frequency mechanical waves with intense oscillations widely used in medicine and cosmetology. By generating a small amount of energy and high-frequency vibrations, these waves create a gentle resonant effect on the skin, improving cell activity, enhancing fat tissue permeability, cleansing sweat glands, and tightening pores.

The ultrasonic facial cleansing scrubber incorporates nanotechnology, generating vibrational oscillations at a frequency of 24,000 Hz per second. This leads to intense molecular friction, accelerating microcirculation and cellular metabolism. The deep high-frequency vibrational massage stimulates microcirculation, promotes fat cell

breakdown, and aids in their removal through metabolic processes.

The spatula-shaped attachment ensures fast and effective pore cleansing by removing dirt, acne, dead skin cells, excess oil, and cosmetic residues.

## Important Safety Measures

- 1.The device was not fully charged at the factory. Charge it before first use.
- 2.Before charging, ensure that the USB charging cable is not damaged or deformed. Do not plug in the device with wet hands.
- 3.Before storing the device, make sure it is turned off. Store it in a dry place, out of reach of children.
- 4.Use antiseptic wipes or a cotton pad soaked in chlorhexidine solution to clean the body and spatula of the ultrasonic scrubber. Do not use corrosive cleaning agents. Clean the device after each use.
- 5.Do not throw the device into fire or store it in high-temperature or humid environments.
- 6.Do not immerse the device in water or rinse it directly under a faucet.
- 7.If the device has been damaged due to a fall, do not use it.
- 8.Do not disassemble or attempt to repair the device yourself. For repairs, please contact the Medica+ service center.

## Indications:

- Dull skin tone
- Enlarged pores
- Oily skin
- Acne and comedones
- Scars and post-acne marks
- Fine wrinkles
- Loss of firmness
- Increased pigmentation

## Contraindications

Although using an ultrasonic scrubber is a relatively safe procedure that can be performed even at home, it is important to remember that it also has several contraindications. The ultrasonic scrubber is not recommended for individuals:

- with couperosis;
- in the acute stage of herpes;
- with a fever;
- with malignant tumors;
- with cardiovascular diseases;

- with arterial hypertension;
- with epilepsy;
- with trigeminal neuralgia;
- during pregnancy (at any stage);
- with increased skin sensitivity;
- with a sunburn or after chemical peeling;
- prone to skin allergic reactions;
- with neurological and dermatological diseases.

## **Usage Recommendations**

Before using the scrubber, remove the protective cap.

Press and hold the power button for 1.5 seconds to turn the scrubber on/off.

Press the button briefly (1 second) to switch between modes and select the desired one (Ultrasonic Cleansing / Lifting Mode / EMS).


# Mode Description

## 1 Mode – Ultrasonic Cleansing

This mode is used for deep facial cleansing, including the removal of blackheads, comedones, pore cleansing, exfoliation of dead skin cells, removal of excess oil and makeup residues, and reduction of pigmentation spots and freckles.

In this mode, hold the scrubber with the sharp side of the spatula against the skin at a 45-degree angle.

Spatula movement – from the periphery to the center.

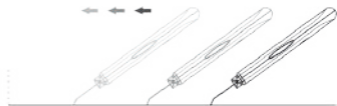
Before the cleansing procedure, it is recommended to apply Medica+ Hydration Gel to the skin for 15 minutes under a film. After this time, remove the film locally and begin the cleansing process. If necessary, additionally moisten the skin with a cotton pad soaked in toner or warm water. After the procedure, apply Medica+ Soothing Mask or Medica+ Pore Control Mask to the skin. For oily and combination skin, the procedure can be performed 1-2 times a week. For dry skin, once a week. The duration of the procedure is individual and depends on the degree of skin contamination. On average, it should not exceed 10-15 minutes.

## Mode 2 – Lifting

This mode is used for lifting and skin tightening, smoothing wrinkles, and enhancing the deeper penetration of skincare products into the skin layers.

In this mode, hold the scrubber with the flat side of the blade against the skin. Move the scrubber slowly from the center to the periphery.

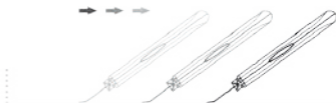The procedure is performed using serums, creams, essences, and facial masks.

This mode can be used daily.

## Mode 3 – EMS

This mode is used for a relaxing facial massage. High-frequency vibrations improve cell activity, boost blood circulation, and accelerate the metabolic process. As a result, the skin becomes firmer, smoother, and softer.

In this mode, hold the scrubber with the flat side of the blade against the skin. Move the scrubber slowly from the center to the periphery.

The procedure is performed using Medica+ Conductive Lifting Gel, serums, creams, essences, and facial masks.

This mode can be used 3-4 times a week.

## Technical Specifications

Product name: Vibroskin 5.0

Battery: 3.7V

Battery capacity: 500mAh

Operating voltage: 3.4V - 4.2V

Rated power: 1.5W

Input voltage: 5V DC

Product dimensions: 194×50×21 mm

Package dimensions: 221×112×40 mm

Net weight: 93g

Frequency: 24 kHz

## Warranty Terms

Dear Customers,

Please keep this manual properly to use the warranty card. The warranty is valid for one year from the date of purchase (or 2-3 years if an extended warranty is purchased) in case of malfunction under normal operating conditions.

The warranty does not cover the following cases:

- If the product is damaged due to improper use or lack of proper maintenance.
- If the product is damaged due to disassembly or modifications.
- If non-original parts from other manufacturers were used with the device.
- Damage caused by force majeure circumstances.

According to the terms, maintenance outside the warranty period is subject to an additional charge.

Parts replaced during the warranty period remain the property of MEDICA+.

## **Package Contents:**

Medica+ Vibroskin 5.0 Ultrasonic Facial Scrubber – 1 unit  
USB Charging Cable – 1 unit

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