# Rotating Electric Toothbrush

ProBrush XO

medica+   
JAPAN TECHNOLOGY

medica+

JAPAN TECHNOLOGY

ENG


**Rotating  
Electric  
Toothbrush  
ProBrush XO**

User Manual

# Content

---

Components .....	01
Technical characteristics .....	02
Before use .....	03
Precautions .....	03
Warning .....	03-04
Charging instructions.....	04-05
Using Guide .....	05-08
Cleaning .....	08
Correct disposal of this product .....	09

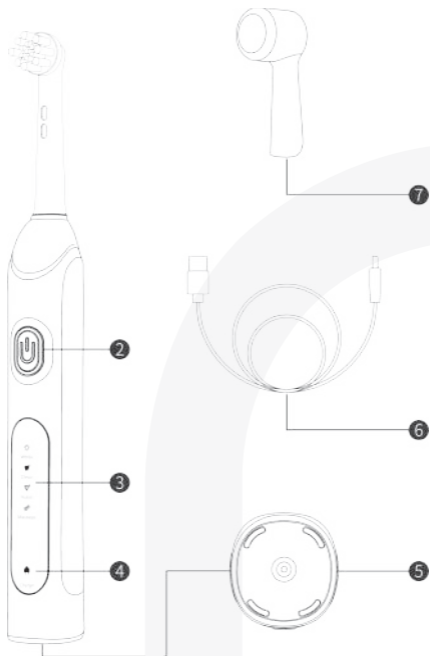1. Brush Head
2. On/Off and Mode selection button
3. Mode indicator
4. Charging indicator
5. Recharge socket
6. USB charging wire
7. Dust cover

Input	DC5V $\approx$ 1.2W
Charging time	About 4 hours
Dimensions	230×26×28mm
Weight	105g

## Before use

Please read this user manual carefully before using and save it for future reference.

## Precautions

Keep the charger away from water. Do not place or store it over or near water contained in a bathtub, washbasin, sink etC. Do not immerse the charger in water or any other liquid. After cleaning, make sure the charger is completely dry before you connect it to the mains.

### Warning

1. This appliance can be used by children and people with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
2. Children should be supervised to ensure that they do not play with the appliance. Cleaning and maintenance shall not be performed by children without supervision.
3. Do not use the charger outdoors or near heated surfaces.
4. Do not clean the brush head, toothbrush handle or charger in the dishwasher.

5. This product is intended to clean teeth, gums and tongue only. Do not use it for any other purpose. Stop using the appliance and consult your doctor for support if you experience any discomfort or pain.
6. Stop using a brush head with crushed or bent bristles. Replace the brush head every 3 months or sooner if signs of wear appear.
7. If the toothpaste you are using contains peroxide, sodium bicarbonate or other bicarbonate (a common ingredient in whitening toothpaste), be sure to wash the brush head thoroughly with soap and water after each use. This prevents possible cracking of the plastic.
8. If you have had oral or gum surgery in the previous 2 months, consult your dentist before using this appliance.

## Charging instructions

1. If the battery is low, charging indicator will flash 10 times to remind user that the battery is low and needs to be charged.
2. Connect the appliance to the DC charging port. Connect the USB cable to any suitable power source. Charging begins with the charging indicator lights up red (a full charge can be completed in about 4 hours). When fully charged, it will auto-off to prevent overcharging and the charging indicator will turn green as reminder.

3. If you use it for the first time or if it is not used for a long time, the indicator may not light up 5 minutes before charging. This the reason for low battery, which is normal.

## Using Guide

1. Insert the brush head

2. Apply the toothpaste to the brush. In order to avoid splashing, guide the brush head to your teeth before switching on the appliance;  
3. Press the on/off button. Press the button again within 3 seconds to choose the desired mode:

### **White Mode:**

Stronger than clean mode. Shocks plaque and dislodges tartar on teeth. Removes 5x plaque in the tooth.

## **Clean Mode:**

Quickly dilutes plaque, and reduces growth of oral bacteria.

(Recommended for everyday brushing)




## **Polish Mode:**

Increase tooth gloss and help remove stains

## **Massage Mode:**

Massage mode to gently stimulate and massage your gums to improve gum health.

4. Guide the brush head slowly from tooth to tooth, spending a few seconds on each tooth surface. Start brushing the outsides, then the insides and finally the chewing surfaces. Brush all four quadrants of your mouth equally (the rotation stops briefly at 30-second intervals to remind). You may also consult your dentist or dental hygienist about the right technique for you.

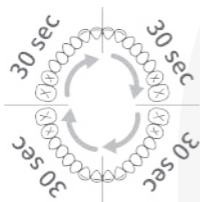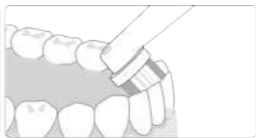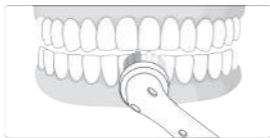The rotation stops briefly at 30-second intervals to remind.

## Note:

In the first days of using any electric toothbrush, your gums may bleed slightly. In general, bleeding should stop after a few days. Should it persist after 2 weeks, please consult your dentist or dental hygienist.

## Cleaning

1. Rinse the brush head with mild waterflow and shake it to remove excess water. You can also take down the brush head, rinse the brush head from the bottom, wash the neck and handle of the toothbrush.

2. Do not clean it in the dishwasher.

## Correct disposal of this product

(Waste Electrical & Electronic Equipment)

This marking shown on the product or its literature, indicates that it should not be disposed with other household wastes at the end of its working life.

To prevent possible harm to the environment or human health from uncontrolled waste disposal, please separate this from other types of wastes and recycle it responsibly to promote the sustainable reuse of material resources.

Household user should contact either the retailer where they purchased this product or their local government office for details of where and how they can take this item for environmentally safe recycling. Business users should contact their supplier and check the terms and conditions of the purchase contract. This product should not be mixed with other commercial wastes for disposal.

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