**ENG / UA / PL**

**SkinLifting PRO 10x**

**medica+**  
JAPAN TECHNOLOGY

medica+

JAPAN TECHNOLOGY

**ENG**

**medica+**  
JAPAN TECHNOLOGY

## **Instruction manual**

**Facial Massager with Microcurrent, RF,  
and Lifting Function**

# **Content**

---


<b>Description .....</b>	<b>01-02</b>
<b>Mode Selection .....</b>	<b>03-06</b>
<b>Operating Instructions .....</b>	<b>07-08</b>
<b>Usage Protocol .....</b>	<b>09-10</b>
<b>Massage Lines .....</b>	<b>11-12</b>
<b>Usage Tips .....</b>	<b>13-15</b>
<b>Contraindications .....</b>	<b>16</b>
<b>Precautions .....</b>	<b>17-18</b>
<b>Warranty Terms .....</b>	<b>19</b>
<b>Technical Specifications .....</b>	<b>20</b>

## Description

The SkinLifting PRO 10x facial massager with microcurrents, RF, and lifting function is a multifunctional cosmetic device for professional skincare at home.

It combines advanced modern technologies to address various skin needs:

- radiofrequency therapy
- electrical muscle stimulation
- dual-action phototherapy.

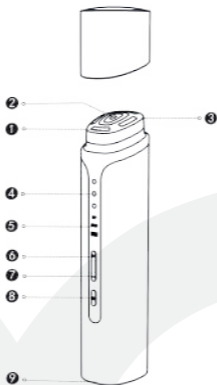1. Electrode plate
2. Temperature sensor
3. Built-in LED light
4. Power/Intensity indicator
5. Mode indicator: Radiofrequency therapy (RF), Electrical Muscle Stimulation (EMS), RF + Vibration
6. Intensity increase button "+"
7. Intensity decrease button "-"
8. Mode switch / Power on & off button
9. Type-C charging port

## Mode Selection

### RF (Radio Frequency) — Radiofrequency Therapy

RF therapy is a procedure that uses radiofrequency waves to heat body tissues. This stimulates the production of collagen and elastin, improving skin elasticity and firmness, reducing wrinkles, and enhancing facial contours. Additionally, RF therapy boosts the absorption and effectiveness of skincare products by promoting deeper penetration of active ingredients into the skin. It also enhances skin radiance, evens out tone and texture, making the face appear fresher and more youthful.

This procedure is recommended for people with dull skin, wrinkles, stretch marks, excess fat deposits, and other signs of skin aging. It is suitable both for preventing age-related changes and for correcting existing skin concerns.

## Electrical Muscle Stimulation (EMS)

This procedure uses electrical impulses to stimulate facial muscles. These impulses cause the muscles to contract and activate, helping to strengthen the facial muscle structure. As a result, facial contours are improved, the jawline is lifted, and wrinkles are reduced. EMS therapy also enhances blood circulation and lymphatic drainage, promoting oxygen and nutrient delivery to skin cells and accelerating the removal of toxins.

It is recommended for individuals experiencing muscle tone loss, sagging facial contours, fine lines, puffiness, and other signs of aging. EMS is suitable both for preventing premature aging and for restoring firmness and skin tone in already aging skin.

## RF (Radio Frequency) + Vibration

A combined treatment that merges the effects of radio frequency waves with mechanical vibration to enhance skin rejuvenation and toning. Radio frequency waves penetrate the deeper layers of the skin, causing controlled heating. As a result, skin elasticity improves, wrinkles are reduced, and facial contours become more defined. Vibration stimulates microcirculation, increases blood flow to the tissues, promotes effective lymphatic drainage, and reduces puffiness. It also relaxes tense muscles, boosts metabolic processes, and enhances the effectiveness of RF energy.

## LED Therapy

- Red light (620–630 nm) stimulates cellular metabolism, activates collagen production, improves blood circulation, and enhances the effectiveness of treatments. It tightens pores, regulates sebum production, increases skin elasticity, and is used in acne treatment.
- Blue light (470–480 nm) works on the skin's surface and has strong anti-inflammatory and antibacterial effects. It eliminates skin dullness and is used in acne and post-acne therapy.

# Operating Instructions

## Charging

1. Before first use, fully charge the device. This ensures proper function, optimal performance, and prolongs battery life.
2. Connect the Type-C charging cable to the corresponding port on the device.
3. When all three intensity/power indicators stay lit — the device is fully charged.

Do not use the device while it is charging. Place the device on its side during charging to prevent it from falling.

## Mode Switching

After turning on the device, briefly press the power/mode button to switch between modes in the following order:  
[RF] → [EMS] → [RF/Vibration].

## Turning the Device On/Off

Press and hold the power/mode button for 2 seconds to turn the device on. By default, the first mode — “RF” — will start. If no action or skin contact occurs within 6 minutes after turning on, the device will automatically power off.

To turn the device off, press and hold the power/mode button for 2 seconds.

## Selecting Intensity

1. By default, the device starts at the first intensity level (1 indicator light).
2. When switching modes, the intensity automatically resets to level 1.
3. Short press the “+” button to increase intensity.
4. Short press the “-” button to decrease intensity.
5. Maximum — 3 levels (3 indicator lights), minimum — 1 level (1 indicator light).

## Battery Level Indicator

The battery level is not displayed during device operation. To check the battery level when the device is off, briefly press the power/mode button:

- 3 indicators — full charge
- 1 indicator — low charge

## Low Battery Reminder

1. When the battery is low, the indicator will flash while the device is on or in standby mode. Please connect the charger as soon as possible.
2. If the device does not turn on and the first indicator blinks 5 times and then turns off — this means the battery is too low and the device needs to be charged.

## Usage Protocol

- **First mode – RF** – recommended to use once every 2 days.
- **Second mode – EMS** – recommended to use once every 2 days.
- **Third mode – RF + Vibration** – recommended to use once every 2 days.

Suggested application sequence for each mode: Left cheek → Right cheek → Forehead

1. Cleanse your facial skin using a foam or gel cleanser.
2. Tone your face with a facial toner.
3. For modes 1 and 2, use Medica+ conductive gel. For mode 3, use a serum/fluid/active ingredient suitable for your skin type.
4. Repeat the procedure for both sides of the face in mirror order. If you feel dryness or resistance while sliding the device — apply more gel or the product you were using.

5. Each mode session lasts approximately 6 minutes, with a reminder every 3 minutes to switch to another zone or side of the face.
6. After completing a session in one mode, if you wish to switch to another mode, turn the device off and then back on again. To change the mode, briefly press the mode switch button.
7. Total treatment time on the face should not exceed 20–25 minutes per day.
8. After the session: Rinse your face with water.
9. Apply a soothing mask or a pore-tightening Medica+ mask for 10–15 minutes, then rinse it off with water.
10. Apply a cream suitable for your skin type.
11. After completing the procedure, gently wipe off any remaining gel from the surface of the device with an antiseptic wipe or a cotton pad moistened with chlorhexidine.

# Massage Lines

Move the device  
along the neck  
from top to  
bottom

Glide from the  
corner of the lips  
toward the ear.

Move upward along  
the nasolabial folds,  
then continue along  
the temples.

Press the device  
gently against the  
skin and slowly glide  
it upward.

- Move the device from the chin to the earlobe with slow, smooth motions.
- Then place the device under the earlobe and move it back and forth in that area.
- Next, move downward along the neck line, avoiding the thyroid gland area.

Glide in a zigzag motion from the eyebrow to the hairline on the forehead, using slow and smooth movements.

## Usage Tips

To treat the left and right cheeks, place the device on the skin and slowly glide it from the center of the face outward. Repeat the motion several times.

In the eye area, reduce the intensity and slow down the movement to avoid skin stretching and irritation.

For the forehead area, place the device on the skin and slowly move it upward — from the bottom of the forehead to the top. During the procedure, do not pause — the device should continuously glide over the skin in slow, smooth motions.

During the first use, be sure to perform a preliminary test to ensure your skin tolerates the procedure well. This helps avoid possible irritation or unwanted reactions.

- Clean the inner part of your forearm and thoroughly dry the skin with a towel.
- Apply a small amount of gel to the cleansed area. Turn on the device and perform the procedure on the skin for 3 minutes. Be sure to test at intensity level 1.
- After completing the test, remove any remaining gel. Then carefully examine the skin for any unusual reactions (redness, itching, rash, etc.).

## Correct Angle and Device Position:

1. Hold the device in your hand and press the electrode head firmly against the skin. Slowly move it across the selected area.
2. Make sure the electrode maintains full contact with the skin throughout the entire session, and that the device is continuously moving — do not leave it in one place.

## Warning:

- To avoid skin issues, always perform a patch test before the first use.
- If you experience itching, redness, rash, or irritation during the test — stop using the device immediately.
- Do not use the device if your skin looks or feels abnormal.

If within 24 hours after the test your skin shows no redness, heat, or discomfort — the device is safe to use.

## Contraindications

- Epilepsy
- Oncological diseases (even in remission — only with doctor's approval)
- Heart conditions or presence of a pacemaker
- Chronic diseases during an exacerbation phase
- Metal implants in the treatment area (e.g., dental implants, plates, etc.)
- Active herpes outbreak
- Damaged skin (wounds, fresh stitches, burns)
- Dermatological conditions (eczema, psoriasis, couperose, rosacea)
- Botox injections — not earlier than 2 weeks after the procedure
- Dermal fillers — generally allowed after 2–4 weeks
- Chemical peels, laser treatments — only after full skin recovery
- Hypersensitive skin or tendency to couperose

## Precautions

1. The device is intended for use only by individuals aged 18 and over. Do not allow children to use it.
2. Do not use the device on animals — this may cause them harm. Use the device strictly as intended and in accordance with the manual, as misuse may result in injury or damage.
3. Individuals with physical or mental impairments should not use the device due to the risk of injury or accidents.
4. Do not use the device if it is malfunctioning — it may pose a health risk. Avoid using the device near flammable substances (alcohol, solvents, gasoline, lacquers, sprays, etc.) or with flammable cosmetics — there is a risk of fire, electric shock, or injury.

5. Do not charge, use, or leave the device near open flames, heaters, or in high-temperature environments — this may lead to overheating, fire, or explosion.
6. Using the device near water, in the bathroom, or while bathing is strictly prohibited — this may result in electric shock, short circuit, or fire.
7. Do not drop or strike the main unit — this may cause damage or injury. If the device is damaged, it must not be used.
8. Avoid contact between the electrodes and metal objects — this is dangerous and may cause malfunction.
9. Always turn off the power during cleaning or when not in use to avoid electric shock or short circuit.
10. Do not disassemble, repair, or modify the device yourself — this may result in fire, electric shock, or injury. In case of malfunction, please contact the Medica+ service center.
11. If you notice any unusual sound, smell, or smoke during charging, immediately disconnect the charging cable and stop using the device.

# Warranty Terms

Dear Customers,  
Please keep this instruction manual in a safe place in order to use your warranty card. The warranty is valid for one year from the date of purchase (or for 2 or 3 years if extended warranty is purchased) in case of malfunction under normal operating conditions.

The warranty does not cover the following cases:

1. If the product is damaged due to improper use or lack of proper maintenance.
  2. If the product is damaged as a result of disassembly or modifications.
  3. If components or accessories not manufactured or approved by our company have been used with the device.
  4. Damage caused by force majeure circumstances.
  5. According to the terms, maintenance services provided outside the warranty are subject to additional charges.
- Any parts replaced during the warranty period remain the property of MEDICA+.

## Technical Specifications

Product Name: SkinLifting PRO 10x Facial  
Massager with Microcurrents, RF, and  
Lifting Function

Size: 160 × 38 × 32 mm

Rated Input: 5V / 1A

Battery Capacity: 2000 mAh

EMS Frequency Modulation: 130 V

Red Light: 650 ± 10 nm

Blue Light: 460 ± 10 nm

RF Frequency (Radio Frequency Energy): 1 ~  
1.3 MHz

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