## Warranty card for service records


The product warranty is provided by the manufacturer for medica+ products from 12-36 months.

## Portable darsonval **Darsoline 7.0**

### User manual
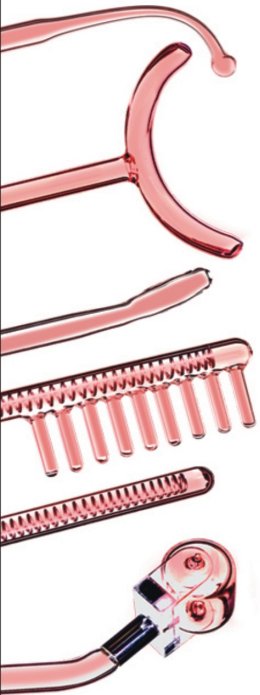
## 1. What is high frequency:

The machine generates a weak high-frequency current through a glass electrode filled with neon or argon using infrared light and ultraviolet light.

A high-frequency electrotherapy device consists of a high-frequency oscillatory circuit, capacitive-resistance and semiconductor devices, equipped with an insertable insulated electrode handle and a glass electrode.

## 2. How does it work?

Inside the handle there are transformers, and on the side of the handle there is an electrical wire for outputting microcurrent.

When you turn on the machine, is there a high frequency high frequency current oscillating circuit? Produces a discharge of a glass electrode and produces infrared or ultraviolet radiation.

When the skin comes into contact with a high frequency current, the vessels will alternately contract and expand.

Electrolysis oxygen sparks from the air and generates ozone, which plays an important role. Has a bactericidal anti-inflammatory effect and improves blood circulation.

The effects of high frequency electrotherapy instruments are as follows: disinfection, acne treatment, acne treatment, reduction and regulation of sebum secretion and metabolism; promotes lymph and blood circulation, improving the ability to regenerate cells. Heat and vibration act on the skin with a smoothing massage.

## 3. Parameter:

Product name:	Portable Darsonval
Model:	Darsline 7.0
Voltage:	110 V / 220 V
Frequency:	50-60 Hz
High frequency vibration:	220 million per second
Power:	10W
Tip material:	Korean ABS plastic
Tube material:	glass
Package size:	38*21*6cm
GW	1 KG

## 4. High frequency can:

1. Sterilize the skin with special ultraviolet / infrared radiation
2. You can select spot machining or large area machining using a different electrode
3. Treat acne
4. Improve hair growth using the high frequency glass hair comb
5. The glass electrode tube can be selected



## 5.7 Electrode specifications:

1. Facial / body massage roller
2. Hair Growth Tube
3. Acupressure massage (large)
4. Bending tube
5. Tongue tube
6. Mushroom tube
7. Fork electrode for large-area massage (body or neck)

Note. Approximate metal head diameter: = 11.55mm (Note: If you intend to use this RF electrode with other RF equipment, please review the above diameter measurements carefully before purchasing to ensure proper compatibility and installation.)
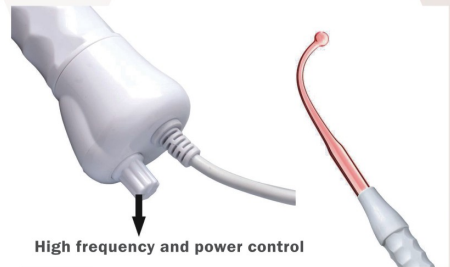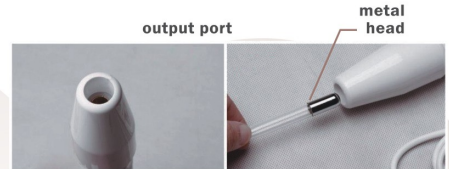
## 6. Procedure before treatment

1. Sterilize the electrode glass with 75% alcohol.
2. Plug in the power supply and insert the stain remover into the appropriate socket.
3. Turn on the power and adjust the density to the desired degree.
4. Take your stain remover and gently hover over the spotted mole until you see the melanin released as waste.
5. Apply disinfectant to the surface.

## 7. Usage

Apply to face or body 3 times a week according to directions.

## Facial / Body Massage Roller

Suitable for larger face and body skin.

- 1) Connect the roller tube to the handle.
- 2) Turn on the machine and microelectricity will appear. pass the.
- 3) Lightly massage your neck or body from the bottom up.
- 4) After use, set the intensity to "OFF". Be sure to turn off the power and remove the glass tube from the handle.

## Fork Electrode for Large Area Massage Tube

Suitable for larger arms, legs, neck and body skin.

- 1) Connect the mushroom tube to the handle.
- 2) Turn on the machine and microelectricity will appear. pass the.
- 3) Lightly massage your neck or body from the bottom up.
- 4) After use, set the intensity to "OFF". Be sure to turn off the power and remove the glass tube from the handle.

## Acupressure massage for face or body (large)

Suitable for treating acne on the face or body.

- 1) Connect the tube to the handle.
- 2) Turn on the machine and microelectricity will appear. pass the.
- 3) Lightly massage your neck or body from the bottom up.
- 4) After use, set the intensity to "OFF". Be sure to turn off the power and remove the glass tube from the handle.

## Mushroom Tube

Suitable for larger areas of the face, cheeks, forehead, neck and body.

- 1) Connect the mushroom tube to the handle.
- 2) Turn on the machine and microelectricity penetrates.
- 3) Lightly massage your face from neck to chin and visor to forehead.
- 4) After use, set the intensity to "OFF" . Be sure to turn off the power and remove the glass tube from the handle.

## Flexible tube

Suitable for parts such as the nose groove. Small parts behind the ear, etc.

- 1) Connect the mushroom tube to the handle.
- 2) Turn on the power and you will see a micro electric current.
- 3) Massage with a glass tube such as "Z" or spiral on the face, the procedure is as follows, before - head-nose-right side of the face-chin-left the face-nose-forehead side.
- 4) Adjust the intensity to "OFF" after use. Be sure to turn off the power and remove the glass tube from the handle

## Language Handset

Is suitable for the skin of the face, for bags under the eyes and for the nose and has the best therapeutic effect.

- 1) Connect the handset to the handle.
- 2) Turn on the machine and the micro-electricity will pass.
- 3) Massage the face from the nose-eyes-right side of the face-chin-left side of the face - eye-nose-forehead.
- 4) Adjust the intensity to "OFF" after use. Be sure to turn off the power and remove the glass tube from the handle.

## Hair Growth Combing Tube

May stimulate the skin surface, reduce folliculitis and aid hair growth

- 1) Connect the comb tube and use it, preferably with some hair restoration product.
- 2) Slowly move the comb tube according to the hair direction.
- 3) Adjust the intensity to "OFF" after use. Be sure to turn off the power and remove the glass tube.

## 8. Main function:

- 1) Helps calm the nerves and has analgesic function.
- 2) Reduces skin irritation, eliminates acne and pimples.
- 3) Accelerates blood circulation and improves metabolism.
- 4) Has anti-aging, anti-wrinkle function and skin tightening.
- 5) Provides ozone sterilization. To calm the nerve with pain relief. To remove blackheads and heal the wound.
- 6) Generation of ozone for sterilization, helping wound healing and restricting pores
- 7) Improved secretion, pH value, causes skin to become neutral.
- 8) Improves blood circulation in the head, stimulates the absorption of nutrients by the skin, treats neurasthenia, headaches, dizziness and insomnia. It also promotes healthy hair growth, prevents and treats hair loss, hair color loss, and more.

## 9. Method of treatment:

### 1). Indirect Method:

- (1). Screw the glass tube into the tip.
- (2). Please ask the customer to hold the glass tube after touching the talcum powder. Turn on the machine and the electricity will flow from your hands to your entire body.
- (3). Cosmetologists massage the client's face from the neck to the chin and from the visor to the forehead. (fourteen). Set the intensity to "0". Make sure the machine is off and pull the glass tube out of the handpiece. This method is suitable for dry natural skin and older skin.

### 2). Direct Method:

- (1). First, lubricate the face with essence and cream, and then the beautician will grab the metal rod to slide the glass tube over the client's face. Turn on the machine.
- (2). Massage with a glass tube type "2" or a spiral on the client's face, the procedure is performed as follows: forehead - nose - right side of the face, chin - left side of the face - + nose-forehead.
- (3). Set the intensity to "0". Make sure the

machine is turned off and remove the glass tube from the handpiece. This method is suitable for oily natural skin.

### **3). Blinking method:**

(1). Cover the patient's eyes with a piece of damp cloth and touch the face with a glass tube from the side of the inflammation. A single treatment should last less than 10 seconds.

(2). It is normal for sparks to appear when the glass tube touches your face.

(3). Adjust the intensity to "0". Make sure the machine is off and pull the glass tube out of the handle. This method is suitable for skin with inflammation and wounds.

### **4). Hair care. Method:**

(1). Use the comb tube with the pilatorium.

(2). In the direction of the hair to move the comb tube.

(3). adjust the intensity to "0". Make sure the machine is off and pull the glass tube out of the handpiece.

This method is designed to stimulate the skin's surface, reduce folliculitis and improve healthy hair growth.

## **10. Caution:**

- 1) Clean and sterilize the stain remover and surface before use.
- 2) After each use, remove the AC adapter and stain remover for cleaning and sterilization.
- 3) When using, do not hit and only for the recommended time to avoid burns on the dermis.
- 4) If the stain does not come off in one use, do not use it again until the scab falls off.
- 5) We recommend that you remove no more than 5 stains at a time.
- 6) After use, the customer is advised not to wet the surface to be treated for 22 hours. Do not eat ginger or anything that is dark in color.
- 7) Do not peel off the scab, let it fall off to avoid infections and burns.
- 8) Do not expose to intense sun exposure for one month.
- 9) For sensitive skin, cover with a piece of veil to avoid over stimulation.
- 10) Do not put on one piece for more than 10 seconds.

11) Pregnant women and people with a pace-maker are not suitable for this treatment.

12) Make sure the glass tube is t inserted into the hand piece.

13) Skin with freckles and splashes is not suitable for this treatment.

## 11. Package List:

Attachment	1pcs
Face / Body Roller	1pcs
Hair Growth Comb for Hair	1pcs
Acupressure, large tube	1pcs
Curl or pinpoint eye Spotting tube	1pcs
Tougue Mushroom Tube	1pcs
Mushroomor Round Tube	1pcs
Fork electrode for large area massage	1pcs
Instructions	1pcs

## 12. Important Safety Notice:

- 1). Please do not disassemble the equipment or try to perform other operations specified in our instructions. All repairs must be carried out by our authorized professional personnel.
- 2). Please do not assemble or use this equipment when you are near water or with wet hands. Please note that do not pour liquid on it.
- 3). Disconnect the plug from the outlet if you have met the following conditions. And stay in touch with a professional service station.
  - a). Liquid has spilled into the equipment.
  - b). smell, smoke, or loud and loud noise.
  - c). The office is broken.
  - d). Equipment fell and broke.
- 4). Please do not place anything on the cables, do not place the equipment in places where someone might step on the cable.
- 5). Unplug the power cord after use for safety reasons.
- 6). Do not insert anything into the equipment slot, or it may cause fire or touch electricity. If something gets caught in the equipment, please

do not remove it yourself, contract the distributor or our company immediately.

7). Please do not place equipment on anything that cannot stand firmly, and avoid tipping over and damaging it.

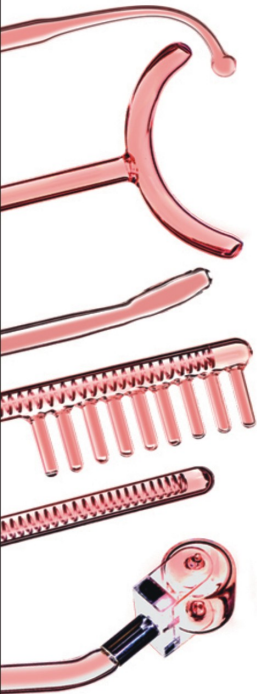
8). Our company reserves the right to change the information in this introduction prior to notice! And leave an explanation to the above.

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