ENG / UA/ PL

PROBRUSH 8.0

medica+
JAPAN TECHNOLOGY

medica+

JAPAN TECHNOLOGY

ENG

Instruction

PROBRUSH 8.0

Contents

Safety Instructions.....	01-02
Description	02-03
Charging	04-05
Replacing the Brush Head	05-06
Brushing Modes	06-07
How to Use	07-08
Cleaning and Storage	08
Warranty Terms	09
Technical Specifications	10

Dear Customer,
Thank you for choosing MEDICA+ products.
Please read this manual carefully before using
the device.


Safety Instructions

1. This gadget is intended for individual use. The ultrasonic toothbrush should not be used by persons with reduced physical, sensory, or mental capabilities unless they are supervised or have been given proper instructions concerning safe usage.
2. Persons who have undergone oral or gum surgery within the past two months are advised to consult a dentist before using the product.
3. Do not use the charger outdoors, in damp environments, or near heat sources or heated surfaces. Avoid any contact of the charger with water.
4. The product is intended solely for hygienic cleaning of teeth, gums, and tongue. Any use of the product for purposes other than its intended function is prohibited. If pain, irritation, or discomfort occurs, discontinue use immediately and consult a doctor or dentist.
5. Do not use brush heads with damaged, deformed, or worn bristles. It is recommended to replace the brush head at least once every 3 months, or earlier if signs of wear appear.

6. When using toothpastes that contain sodium bicarbonate, hydrogen peroxide, or other abrasive ingredients (common in whitening toothpastes), thoroughly rinse the brush head with running water after each use.
7. The manufacturer is not liable for any damage to the device or harm to the user resulting from improper use or failure to follow the instructions provided in this manual.

Description

The MEDICA+ PROBRUSH 8.0 ultrasonic toothbrush operates at a frequency of up to 40,000 vibrations per minute, providing significantly more effective plaque removal compared to manual brushing. The toothbrush features 5 cleaning modes and is equipped with a built-in 2-minute timer with 30-second intervals to help maintain the dentist-recommended brushing duration.

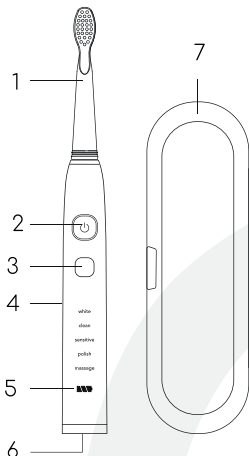- 1 — Toothbrush head
- 2 — Power button
- 3 — Mode selection button
- 4 — Toothbrush handle
- 5 — Charging indicator
- 6 — Charging port
- 7 — Travel case

A tickling or tingling sensation may occur during initial use due to the sonic vibrations. This is a normal reaction and typically disappears as you become accustomed to the device.

Charging

Before first use, charge the ultrasonic toothbrush for 8 hours or until the battery is fully charged. If the red charging indicator is flashing or the device does not turn on, this indicates a low battery level and the need to connect the device to the charger.

If the battery has not been charged for a long time, the indicator may not light up during the first 5 minutes of charging. This is normal.

While charging, the indicator glows red, indicating a low or medium battery level. Once fully charged, the device automatically stops charging to prevent overcharging (see the table below).


Indicator Status	Battery Charge Level
Indicator is off, device not working	< 10 %
Indicator flashes red	10 % – 30 %
1/3 of the indicator flashes green	30 % – 50 %
2/3 of the indicator flashes green	50 % – 80 %
Full indicator flashes green	80 % – 100 %
Indicator turns off	100% (fully charged)

The toothbrush does not operate while charging.

Replacing the Brush Head


Installation

Insert the toothbrush head into the handle of the device. Align mark "A" with mark "B" as shown in the diagram. Rotate the brush head until mark "A" is fully aligned with mark "T".

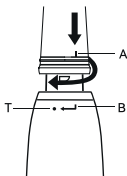Removal

Rotate the toothbrush head until mark "A" aligns with mark "B" as shown in the diagram, then gently remove the head by pulling it along the axis of the handle.

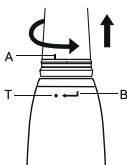Do not use brush heads with damaged or deformed bristles, as they may break during use. It is recommended to replace the brush head every 3 months.

Brushing Modes

When the toothbrush is turned on for the first time, it starts in White mode. Pressing the "mode" button again quickly switches to the next mode.

The brush has five operating modes:

White – 40,000 vibrations/min

Intensive cleaning and removal of surface stains, ideal for individuals without increased tooth sensitivity.

Clean – 36,000 vibrations/min

Gentle daily cleaning, suitable for sensitive teeth and gums.

Sensitive – 34,000 vibrations/min

Very delicate mode for people with highly sensitive teeth and gums, providing soft cleaning without irritation.

Polish – 30,000 vibrations/min

Mode for gentle polishing of the tooth surface, helping to maintain enamel smoothness and a clean feel without putting stress on teeth and gums.

Massage – 32,000 vibrations/min



Mode for gentle gum massage and stimulation, promoting healthy circulation and overall oral care.

To turn off the toothbrush, press and hold the power button for 3 seconds after use. The brush remembers the last used mode and will start with it the next time.

How to Use

Dampen the brush head, apply a small amount of toothpaste, and turn on the device by pressing the power button after placing the brush inside your mouth. Use the “mode” button to select the desired cleaning mode.

While brushing, position the brush head at a 45° angle to the surface of the teeth and gum line. Gently and slowly move the brush along the dental arch without applying additional pressure. Ensure the bristles are in contact with both the teeth and gums at the same time, focusing on the area between them for maximum cleaning effectiveness.

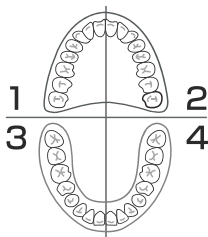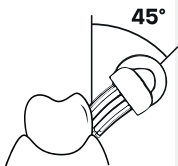It is recommended to mentally divide the oral cavity into four zones and clean each of them for 30 seconds. The intelligent timer automatically tracks the total brushing time and ensures compliance with the dentist-recommended two-minute oral hygiene routine.

Cleaning and Storage

- After each use, thoroughly rinse the brush head and bristles under running water. Wipe the brush handle with a towel or cloth.
- For optimal performance, replace brush heads every 3 months or sooner if signs of wear appear—only use original Medica+ accessories. Do not immerse the charging device in water or expose it to running water.
- Clean the toothbrush and store it in a cool, dry place, away from direct sunlight and out of reach of children.

Warranty Terms

Dear Customers,

Please keep this manual in a safe place in order to use your warranty certificate. The product is covered by a one-year warranty from the date of purchase (or 2 to 3 years if extended warranty was purchased), in the event of a malfunction under normal operating conditions.

The warranty does not cover the following cases:

- If the product is damaged due to improper use or lack of proper maintenance.
- If the device has been disassembled or modified.
- If non-original parts not supplied by our company have been used.
- Damage caused by force majeure circumstances.

According to the policy, service outside the warranty period is subject to additional charges.

Any parts replaced during the warranty period remain the property of MEDICA+.

Technical Specifications

- Power: 1 W
- Input Voltage: 5 V
- Input Power: 1.5 W
- Noise Level: ≤ 65 dB
- Charging Time: ≤ 8 hours
- Vibration Frequency: 30,000–40,000 vibrations/min (depending on selected mode)
- Bristles: DuPont™
- Bristle Hardness: Soft / Medium
- Battery Life: Up to 200 days on a single charge
- Battery Capacity: 1500 mAh
- Weight: 110 g
- Water Resistance: IPX7
- Dimensions (H × W × D): 29 × 27 × 23.6 mm
- Materials: ABS / TPE / Nylon

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