ENG / UA/PL

FootMass 8.0

medica+ 
JAPAN TECHNOLOGY

ENG

Instruction

FootMass 8.0

Contents

Safety Instructions	01-02
Indications for Use	02
Product Description	03-04
Screen Overview	04
Button Overview.....	05-08
Wireless Remote Control Overview	09
Device Usage	10-11
Device Cleaning	12
Troubleshooting	14
Warranty Policy	15
Technical Specifications	16

Dear Customer,

Thank you for choosing MEDICA+ products.

Please read this manual carefully before using the device.

Safety Instructions

1. Monitor the usage time — do not exceed 30 minutes per session. Exceeding this limit may cause muscle or nerve discomfort.
2. Do not use the product on skin with tumors, inflammation, or open wounds.
3. Do not use if you have acute painful conditions such as osteoporosis, spinal stenosis, fractures, muscle tears, or muscle detachment.
4. If you experience any discomfort during use, stop using the device immediately.
5. Not recommended for children or individuals with limited self-care abilities or sensory impairments.
6. Persons with physical or mental disorders should use the device with caution and only under supervision.
7. Keep away from open flames — this may cause an accident or fire hazard.
8. Do not use the device in humid environments (bathroom, kitchen, near water).

9. Do not allow the device to fall into a bathtub, sink, or other containers filled with water. Do not immerse the device in water.

Indications for Use

1. Poor blood circulation in the legs due to prolonged sitting, improper muscle load, or tight footwear.
2. Recovery after physical activity or sports training.
3. Sleep disturbances related to stress or tension in the legs.
4. Need for metabolic stimulation and detoxification to support immune function.
5. Chronic pain or discomfort in the feet, including plantar fasciitis and neuropathy.
6. Symptoms of depression or low mood that may be alleviated by stimulating foot reflex points.
7. Swelling of the feet and ankles associated with fluid retention.

Product Description

The MEDICA+ FootMass 8.0 foot massager combines rolling and air compression massage with a heating function. It effectively helps relieve various types of foot pain, including neuropathic numbness, chronic discomfort caused by plantar fasciitis, and reduces tension while promoting deep relaxation. Additionally, it improves blood circulation and supports healthy lymphatic drainage.

The best time to enjoy a foot massage is after sports activities, exercise sessions, or long hours in high heels during a busy workday or walk. With this powerful device, you can experience both gentle and deep massage, delivering a profound sense of comfort and relaxation.

Removable
(washable) cover


Control panel

LCD display

Leg cuffs

Screen Overview

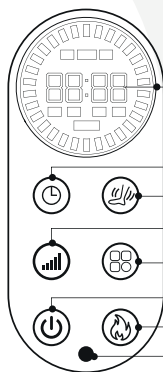LCD display

Timer

Foot massage mode

Intensity

Mode selection button

Power on/off button

Heating

Remote control

Button Overview

Power Button

Controls the start and shutdown of the massager.

 **Heat.** Used to control the heating function. A short press turns on the heat, and the word “HEAT” appears on the display. Two heating levels are available and can be toggled with short presses of the same button. To turn off the heat, press the button again — the “HEAT” indicator will disappear from the screen.

 **Timer.** Allows you to set the massage duration for each session. Press briefly to choose one of six time intervals: 5, 10, 15, 20, 25, or 30 minutes. The selected time will be shown on the display. Once the countdown finishes, the device will automatically turn off.

 **Foot Massage Mode.** Activates the air-compression massage for the legs. By default, the mode is off. A short press turns it on, and the display will show “90”.

There are three intensity levels — F1, F2, and F3, which can be cycled through with short presses of the same button. To turn off foot massage mode, press the button again — the function will deactivate, and the display will clear.

Note: Before using the foot massage function, connect the air tube to the socket located at the back of the device, as shown in the image below. Also, attach the compression cuff to enable air-compression massage. This enhances the functionality of the device and increases the effectiveness of your massage sessions.

⊞ Mode Selection Button

Five massage modes are available:

Mode 1: Airbag massage + roller rotation for 60 seconds at low frequency, with periodic foot hold. The roller alternates direction (forward/backward) every 60 seconds.



Mode 2: Airbag massage + roller rotation for 30 seconds at medium frequency, with periodic foot hold. The roller changes direction every 30 seconds.

Mode 3: Airbag massage + roller rotation for 10 seconds at medium frequency, with periodic foot hold. The roller reverses direction every 10 seconds.

Mode 4: Airbag massage off. The roller rotates for 60 seconds at low frequency, reversing direction every 60 seconds.

Mode 5: Airbag massage on, roller movement off.

In Mode 4, only the rollers are active, the airbag function is disabled, and "AIR" will not appear on the screen.

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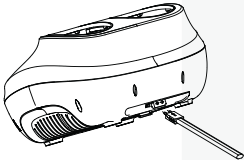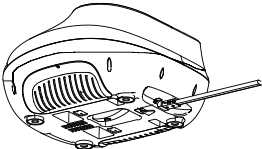 **Intensity.** Select the air pressure strength of the airbag massage. There are 9 adjustable levels, from L1 to L9, with gradually increasing pressure.


Default setting: L5 (automatically set when the device is turned on). Massage intensity can be adjusted to your personal preference:

- L1-L3: Gentle, delicate massage — recommended for elderly users.
- L3-L6: Medium intensity — ideal for most users.
- L7-L9: Deep tissue massage — for those who prefer strong pressure and stimulation of acupuncture points.

Warning:

Do not step on the device during or after a massage session — this may cause damage to the device or result in a fall due to its unstable surface.

Wireless Remote Control Overview

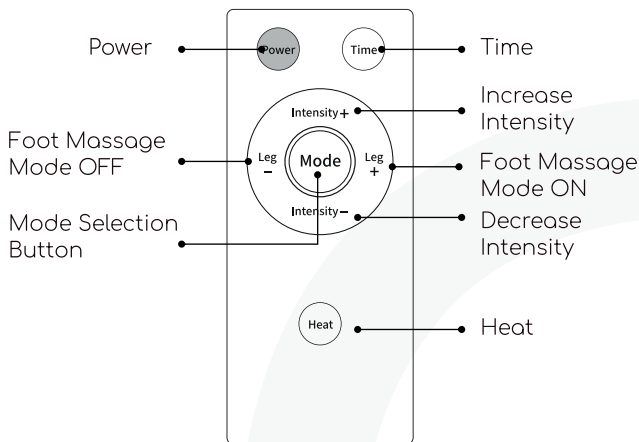If bending down during use is uncomfortable, you can operate the device using the wireless remote control. The remote control buttons fully correspond to those on the LCD control panel.

- The remote works effectively within a range of up to 1.5 meters.
- To use it, you need to install one AAA battery (not included).

Device Usage

1. Before use, ensure the device is not damaged. Then plug the power cord into an outlet, making sure the massager is placed on a stable surface.
2. When a blue indicator light appears on the  display, it means the device is powered on.
3. After placing your feet into the massager, briefly press the Power button  to start the device. By default, Mode 1 is activated, with a 15-minute massage duration and airbag intensity level L5.
4. During use, you can Press the Mode button  to select your preferred massage program. Press the Intensity button  to adjust the air compression level. Press the Foot Massage button  to activate or deactivate the foot air massage. Press the Heat button  to toggle the heating function. Press the Time button  to set the massage duration.
5. When the set time elapses, the massager will automatically shut off. If you wish to stop the massage earlier, press the Power button  to turn off the device manually.

On the first use, the heat function is enabled by default, while vibration is turned off. If you experience any discomfort due to pressure or massage mode, discontinue use immediately.

- The recommended massage duration is no more than 30 minutes per day.
- When adjusting the timer, the countdown updates automatically. If you change the massage mode during operation, the timer continues based on the previously set duration.

Cleaning and Storage of the Device

Cleaning

After each use, make sure to turn off the device and unplug it from the power outlet. If the device becomes dirty, wipe its surface with a damp, soft cloth. Do not use any chemical agents such as gasoline, alcohol, or harsh cleaners. After cleaning, dry the device with a clean cloth and let it air-dry naturally. Do not use a hair dryer or any artificial heat sources for drying.

Removing and Washing the Fabric Cover

Open the zipper.

Remove the cover.

Reinstalling the Cover

Lay out the fabric cover.

Close the zipper.

Storage

1. Unplug the device from the power outlet after each use.
2. Store the massager in a cool, dry place away from direct sunlight.
3. Do not place heavy objects on top of the device to avoid damage.
4. Avoid placing sharp objects on the device to prevent scratches or surface damage.
5. Keep the device out of reach of children to prevent accidental collisions or misuse as a toy.
6. Do not place items containing liquids on top of the packaging box — this may cause leakage into the device and damage internal electronic components.

Troubleshooting

Issue	Possible Causes	Solution
Device does not work	The power plug is not connected to the outlet.	Check whether the local voltage matches the voltage indicated in the specifications and plug the device into a power outlet.
	The device has just completed a 15-minute massage session and has entered standby mode.	If you wish to start another session, please press the "Power" button again.
	The device was used for 30 minutes, which triggered the overheating protection or automatic power-off.	Unplug the device from the outlet and allow it to cool for 60 minutes before using it again.
Device overheating	The heating function is turned on.	Turn off the heating function.
	The device has been used for 30 minutes.	The device has been used for 30 minutes. Disconnect the power supply and allow the device to cool down. Always unplug the device after each use.
Sound from the Massager	The sound is caused by the mechanical structure of the device.	This is a normal operating noise and does not indicate a malfunction.

Warranty Policy

Dear customers,

Please retain this manual properly, as it serves as your warranty certificate. The product is covered by a one-year warranty from the date of purchase (or 2-3 years if extended warranty was purchased), provided it is used under normal operating conditions.

The warranty does not cover the following cases:

1. Damage caused by improper use or lack of proper maintenance.
 2. Damage resulting from disassembly or unauthorized modifications.
 3. Use of parts or accessories not supplied or approved by our company.
 4. Damage caused by force majeure circumstances.
- According to the policy, any maintenance or service outside the warranty scope will be subject to additional charges.

Any parts replaced during the warranty period remain the property of medica+

Technical Specifications

Power: 60 W

- Product Dimensions: 47.5 × 27 × 43 cm
- Material: ABS + Fabric
- Voltage: 100-240 V
- Massage Modes: 5 modes
- Intensity Levels: 9 levels
- Heating Levels: 2 levels
- Heating Temperature: 40°C to 50°C
- Rated Voltage: 100-240 V
- Timer Settings: 5 / 10 / 15 / 20 / 25 / 30 minutes
- Suitable Foot Size: EU 35 to 48
- Cord Length: 2 meters
- Cuff Height: 26 cm
- Cuff Diameter: 18 cm

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